
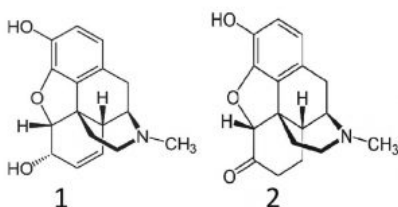


Analysis of Isobaric Opioids by LCMS – AppNote

Separation of Isobaric Drugs

Morphine and Hydromorphone are Isobaric Drugs (the same m/z value) that can be separated using the Cogent Bidentate C18 2.0™ Column. Multiple samples (n=5) were prepared and analyzed. The results showed great reproducibility (RSDs < 3%). The detection limits were 50 to 500 times lower than that required by immunoassays.

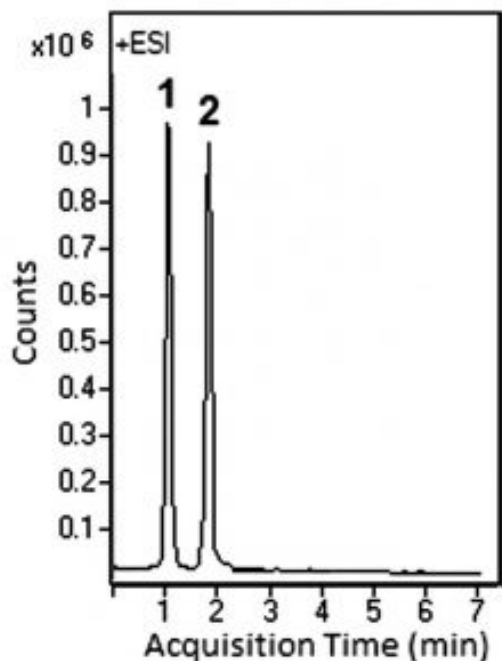

Peaks:

1. Morphine 286.1438 m/z [M+H]⁺
2. Hydromorphone 286.1438 m/z [M+H]⁺

Method Conditions

Column: Cogent Bidentate C18 2.0™, 2.2μm, 120Å

Catalog No.: 40218-05P-2

Dimensions: 2.1 x 50 mm

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: 50% Acetonitrile / 50% Methanol / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
----------------	----

0	5
4	50
5	90
6	90
7	5

Post Time: 3 minutes

Injection vol.: 1µL

Flow rate: 0.4mL/minute

Detection: ESI - POS - Agilent 6210 MSD TOF Mass Spectrometer

Sample: The drugs were spiked into urine at a level of 50 ng/mL. Extraction method: Spiked urine sample was loaded into SPE cartridge I (Clean Screen Xcel™, UCT Bristol, PA, USA) and eluted with 0.78 mL of Acetonitrile, 200 µL of 2-Propanol, and 20 µL of Ammonia. After the elution, the sample was dried under N2 gas and dissolved in 100 µL of 50% Methanol / 50% DI Water / 0.1% Formic Acid. Before injection, the sample was filtered through a 0.45 µm Nylon Syringe Filter (MicroSolv Tech Corp.).

Note: Opioids have been used as pain relieving drugs but are also abused as drugs producing feelings of euphoria. They are highly addictive. Nearly three out of four drug overdoses involve pain killers, so the screening of this class of drugs in urine or plasma samples is very important.



Attachment

No 312 Morphine and Hydromorphone in Urine pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Date: 05-20-2024