
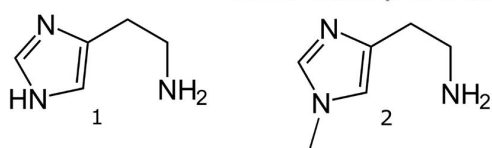


## Histamine and methylhistamine analyzed with LCMS – AppNote

### No derivatization required for detection or retention with this method

Various assay methods for histamine (HA) and / or its metabolite, methylhistamine (MHA), in biological samples, have been developed. However, most of them require post-column (*for detection purposes*) or pre-column (*to achieve retention*) derivatization.

The method in this application note provides separation and detection of these two compounds yet doesn't require any derivatization. It is able to solve the inherently difficult problem of analyzing two biogenic amines with close physicochemical properties. A successful validation of the assay was indicated by the high linearity of calibration curves and the low inter- and intraday variation coefficients.

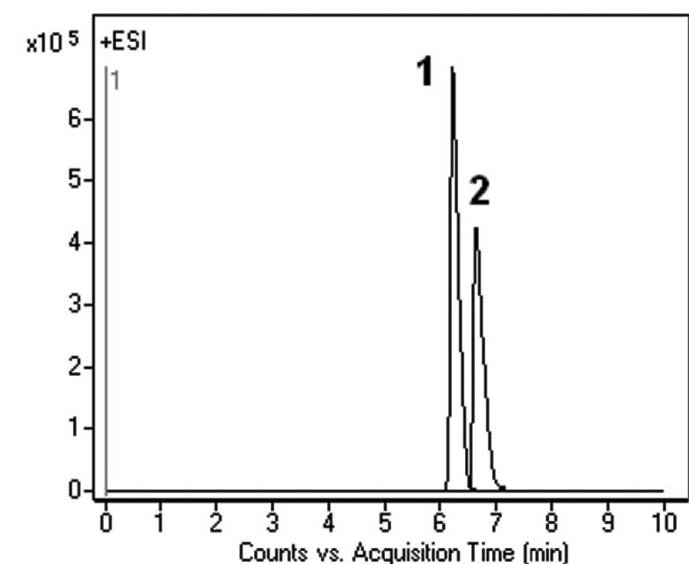

### Peaks:

1. Histamine 112.0869 m/z [M+H]<sup>+</sup>
2. Methylhistamine 126.1026 m/z [M+H]<sup>+</sup>

### Method Conditions

**Column:** Cogent Diamond Hydride™, 4μm, 100Å

**Catalog No.:** 70000-15P-2

**Dimensions:** 2.1 x 150mm

**Mobile Phase:**

A: DI water with 0.1% formic acid (v/v)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

B: Acetonitrile with 0.1% formic acid (v/v)

**Gradient:**

| Time (minutes) | %B |
|----------------|----|
| 0              | 70 |
| 2              | 65 |
| 6              | 10 |
| 8              | 10 |
| 9              | 70 |

**Post Time:** 2 minutes

**Injection vol.:** 1µL

**Flow rate:** 0.4mL / minute

**Detection:** ESI - POS - Agilent 6210 MSD TOF Mass Spectrometer

**Sample Preparation:** 1mg / mL of each compound was dissolved in DI water and filtered through a 0.45µm syringe filter (MicroSolv Tech Corp). Sample for injection was diluted 1:100 with 50:50 solvent A / solvent B mixture.

**t<sub>0</sub>:** 0.9 minutes



**Attachment**

**No 197 Histamine and methylhistamine analyzed with LCMS** [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)