
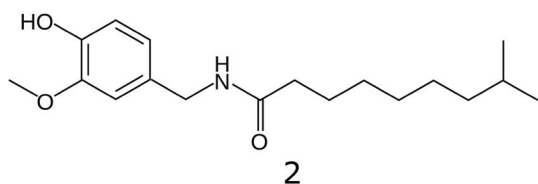
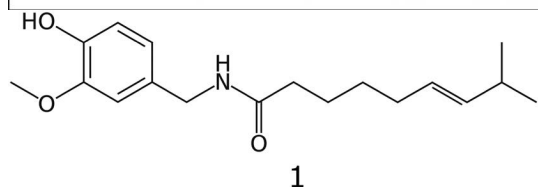


Capsaicin in Topical Solution Analyzed with HPLC – AppNote

Separation of Capsaicinoids in an OTC Formulation

A formulation of Capsaicins is used as a topically applied solution for Arthritis pain relief. Analyzing Capsaicin and a similar compound, Dihydrocapsaicin is required. The Method illustrates the ability to resolve structurally similar Capsaicinoids from a pharmaceutical formulation. Peak Shapes and Efficiencies are excellent for both Peaks.

Furthermore, this AppNote presents data from two Column lots with the overlay of two runs in the figure below to exhibit Robustness and Precision of the Method.

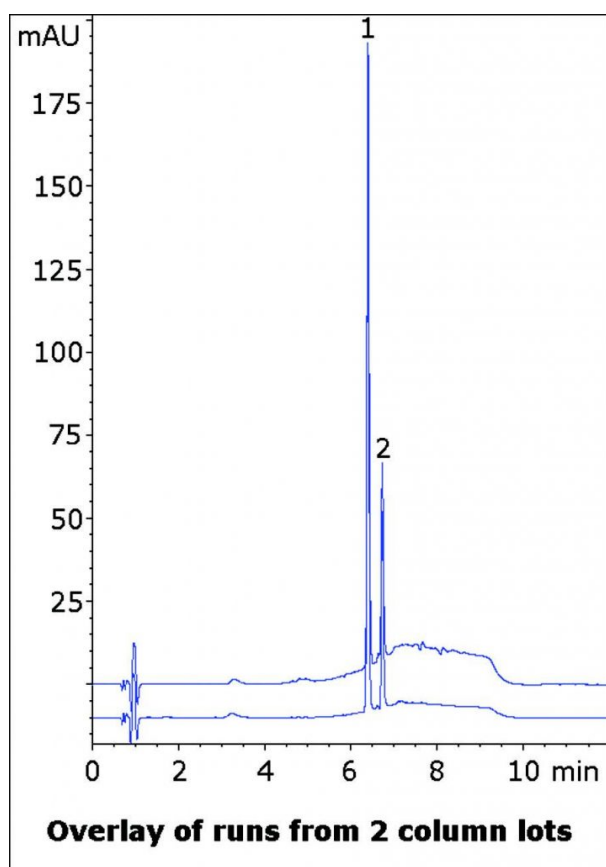

Peaks:

1. Capsaicin
2. Dihydrocapsaicin

Method Conditions

Column: Cogent Bidentate C18™, 4µm, 100Å

Catalog No.: 40018-75P

Dimensions: 4.6 x 75mm

Mobile Phase:

A: DI Water with 0.1% Formic Acid (v/v)

B: Acetonitrile with 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	20
1	20
3	80
8	80
9	20

Post Time: 3 minutes

Injection vol.: 1µL

Flow rate: 1.0mL / minute

Detection: UV @ 280nm

Sample Preparation: TopCare® Solution, a commercially available Pharmaceutical Product, containing 0.15% Capsaicin was filtered with a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.). Peak identities were confirmed using individual USP Reference Standards.

t₀: 0.9 minutes

Note: Capsaicin is the active ingredient in Chili Peppers and has many uses. It is commonly used in foods to add spiciness, but is also used in pepper spray. It is a highly potent irritant in mammals and pure Capsaicin reference standards should be handled with care.



Attachment

No 222 Capsaicin in Topical Solution Analyzed with HPLC pdf 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Date: 09-26-2024