
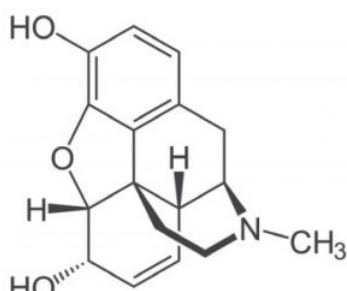


Morphine Sulfate Tablet Analyzed by HPLC- AppNote

Excellent Peak Shape for Tertiary Amine

Morphine can be a very problematic compound for analysis by HPLC. It has a tertiary amine and has a tendency to tail under some conditions. USP assay methodology uses phosphates, alkylsulfonates, or other LC-MS incompatible agents in the Mobile Phase.

This Method uses only Formic Acid and can be readily adapted for LCMS. As such, analysis of Morphine in blood samples can be achieved with excellent retention and peak shape.

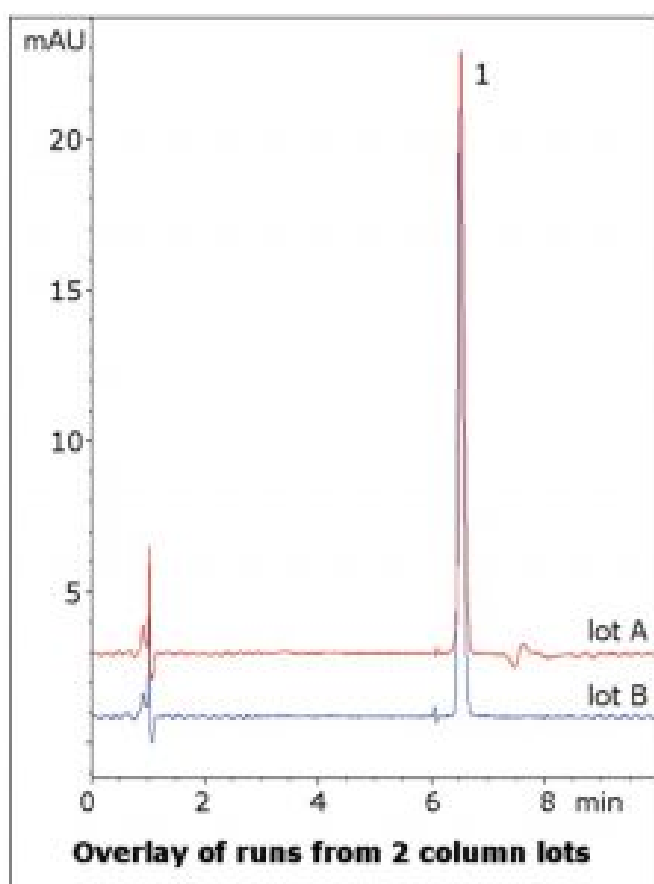

Morphine Sulfate

Peak:

Method Conditions

Column: Cogent Diamond Hydride™, 4µm, 100Å

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Catalog No.: 70000-7.5P

Dimensions: 4.6 x 75mm

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	90
1	90
6	40
7	90

Post Time: 3 minutes

Injection vol.: 1µL

Flow rate: 1.0mL / minute

Detection: UV @ 284nm

Sample Preparation: 15mg strength Morphine Sulfate tablet was ground and weighed in a 25mL volumetric flask. A portion of 50:50 Solvent A / Solvent B diluent was added and the flask was sonicated 10 minutes. It was then diluted to mark and filtered with a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.).

t₀: 0.9 minutes

Note: Morphine is a highly potent opiate analgesic widely used in clinical applications to treat severe pain. However, tolerance and addiction develop rapidly with its use so it has potential for abuse as well. It is named after Morpheus, the Greek god of dreams.



Attachment

No 230 Morphine Sulfate Tablet Analyzed by HPLC pdf 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com