
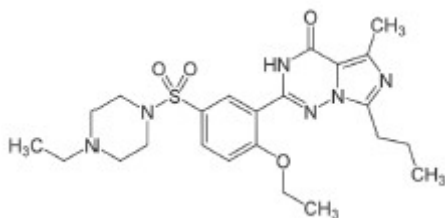


## Vardenafil (Levitra®) Analyzed with HPLC – AppNote

### LC-MS Compatible Assay Method

Vardenafil in a tablet formulation can be readily assayed with this LC-MS compatible gradient method. The peak tailing factor was close to unity. The compound has several amine groups which can produce tailing with ordinary HPLC Columns that have a number of surface Silanols.

The MS-compatible Mobile Phase means that the method can be adapted to more complex samples such as plasma. Five runs are shown in the figure to illustrate the repeatability of the data.

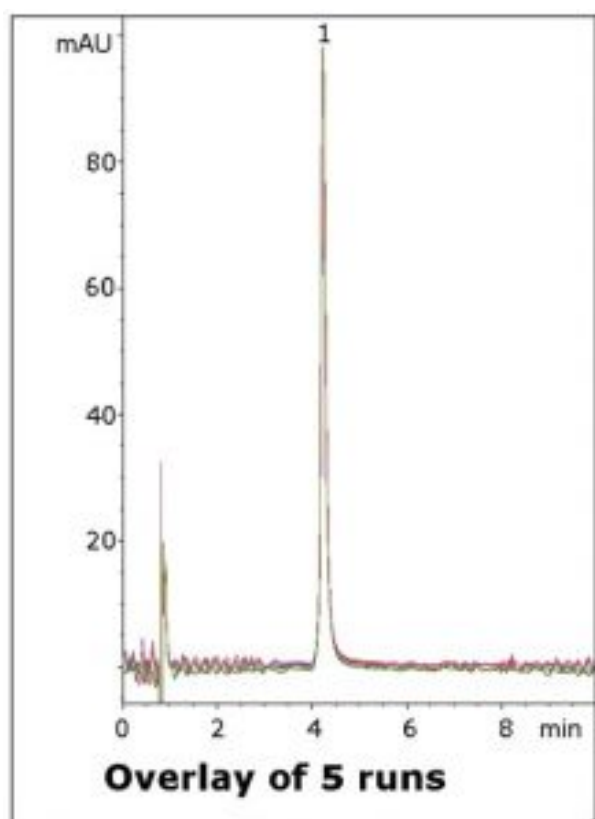

**Peak:**

Vardenafil

### Method Conditions

**Column:** Cogent UDA™, 4µm, 100Å

**Catalog No.:** 40031-7.5P

**Dimensions:** 4.6 x 75 mm

**Mobile Phase:**

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

**Gradient:**

Time (minutes)	%B
0	90
1	90
6	40
7	90

**Post Time:** 3 minutes

**Injection vol.:** 2 $\mu$ L

**Flow rate:** 1.0 mL / minute

**Detection:** UV @ 210nm

**Sample Preparation:** 20mg strength Levitra® tablet was ground and added to a 25mL volumetric flask. A portion of 50 / 50 Solvent A / Solvent B diluent was added and the flask was sonicated 10 minutes. Then it was diluted to mark and mixed. A portion was filtered through a 0.45 $\mu$ m Nylon Syringe Filter (MicroSolv Tech Corp.).

**t<sub>0</sub>:** 0.7 minutes

**Note:** Vardenafil is used to treat Erectile Dysfunction and acts by inhibition of Phosphodiesterase Type 5 (PDE5). It is sold under trade names Levitra® and Staxyn®.

**Attachment**

No 266 Vardenafil by HPLC pdf 0.4 Mb [Download File](#)

Printed from the Chrom Resource Center

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

Date: 05-19-2024