
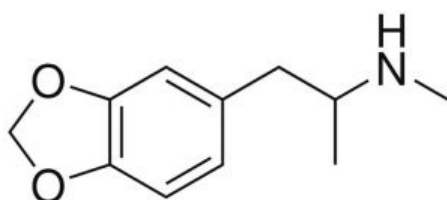


Analysis of MDMA in Plasma Samples with LCMS – AppNote

Methylenedioxymethamphetamine Analyzed with MS

Click [HERE](#) for Column Ordering Information.

Under the described conditions, MDMA was retained and eluted as a Symmetrical Peak. The Sensitivity of the Method is very good and comparable to that reported with GCMS Detection [1]. Matrix effects were of minor extent and reproducible and hence should not compromise Quantification. The Method can be used for Forensic Research and Clinical Analysis.

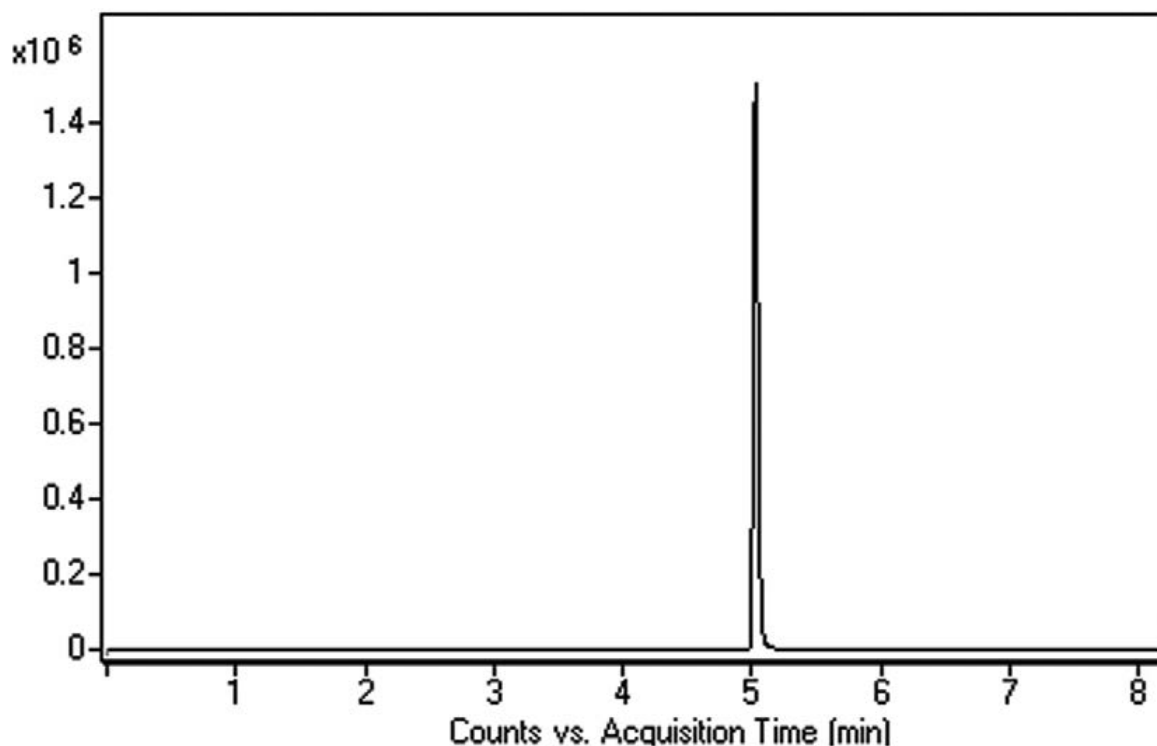

Peak:

(±)-3,4-Methylenedioxymethamphetamine, m/z 194.1176 [M+H]⁺

Method Conditions

Column: Cogent Phenyl Hydride™, 4 μ m, 100Å

Catalog No.: 69020-05P-2

Dimensions: 2.1 x 50mm

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	10
1	10
6	90
7	10

Post Time: 3 minutes

Flow rate: 0.4mL / minute

Injection vol.: 1µL

Sample Preparation: 50 µl of Acetonitrile was mixed with 50µl of plasma for protein precipitation. The samples were centrifuged (16000×g for 15 minutes), and the supernatant was filtered through a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.) and transferred to autosampler vials for injection.

Detection: ESI - POS - Agilent 6210 MSD TOF Mass Spectrometer

t₀: 0.9 minutes

Note: The Amphetamine derivative 3,4-methylenedioxyamphetamine (MDMA), known also as Molly or Ecstasy, is often used or abused as a recreational drug. Because of a reported high inter-individual difference of its toxicity, sensitive analytical methods are needed. A urine test is a standard method to investigate drug abuse but the method has a very low diagnostic sensitivity and makes testing in plasma much more suitable.

Reference:

[1]. R. Kikura, Y. Nakahara, T. Mieczkowski, F. Tagliaro, *Forensic Sci. Int.* 84 (1997) 165-177.

**Attachment**

No 263 Analysis of MDMA in Plasma Samples with LCMS pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com