
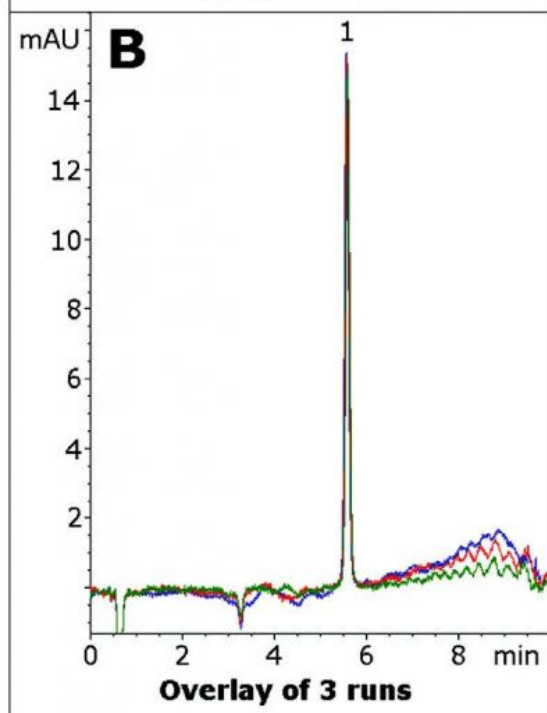


## Vancomycin HPLC method transferred – AppNote

### Benefits of using 2.2 $\mu\text{m}$ column for Vancomycin analysis

The highly polar antibiotic Vancomycin can be readily analyzed with either the standard 4 $\mu\text{m}$  Cogent Diamond Hydride or the Cogent Diamond Hydride 2.0 column. Notably higher efficiency is obtained on the smaller particle size column.

This Method is easy to perform and is LCMS compatible. Three runs were performed on each column in order to demonstrate consistency and show the benefits of using the smaller particle column.

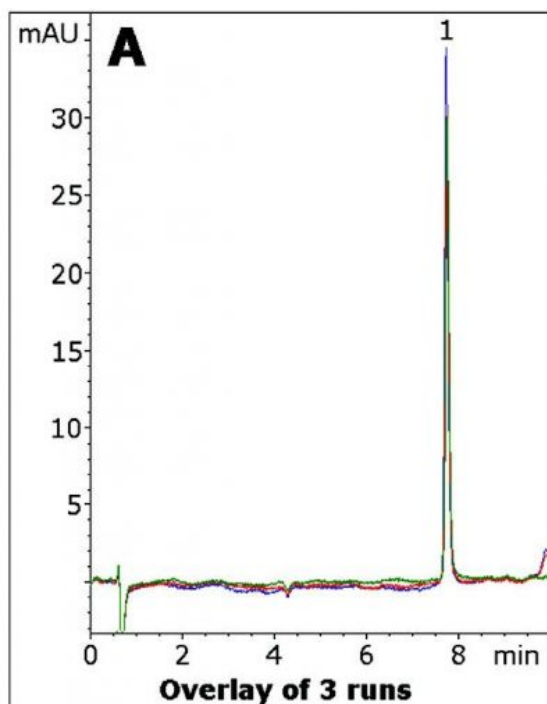

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply


**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

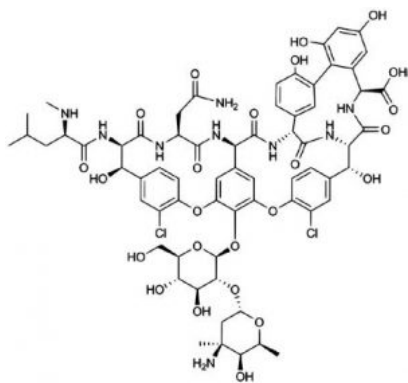

## Peak:

Vancomycin

## Method Conditions

### Columns:

Fig. A: Cogent Diamond Hydride 2.0™, 2.2µm, 120Å

Fig. B: Cogent Diamond Hydride™, 4µm, 100Å

### Catalog Nos.:

Fig. A: [70200-05P-2](#);

Fig. B: [70000-05P-3](#)

### Dimensions:

Fig. A: 2.1 x 50 mm

Fig. B: 3.0 x 50 mm

### Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Time (minutes)	%B
0	70
0.4	70
7.4	10
8.4	70

**Post time:** 5 minutes

**Injection vol.:** 0.2 µL

### Flow rate:

Fig. A: 0.29 mL / minutes

Fig. B: 0.50 mL / minutes

**Detection:** UV @ 210 nm

### Sample:

**Stock Solution:** 1 mg / mL Vancomycin HCl in 50 / 50 solvent A / solvent B diluent. The solution

was filtered through a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.).

**Working Solution:** Stock solution was diluted 1:100 with 50/50 solvent A/ solvent B mixture.

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

**t<sub>0</sub>**: 0.4 minutes

**Note:** Vancomycin is a glycosylated non-ribosomal Peptide antibiotic used to treat Colitis. Vancomycin is often used as a drug of last resort when other antibiotics are rendered ineffective due to developed resistance of bacteria. It is a natural product isolated from *Amycolatopsis Orientalis*.



## Attachment

**Vancomycin Method Transfer pdf** [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)