
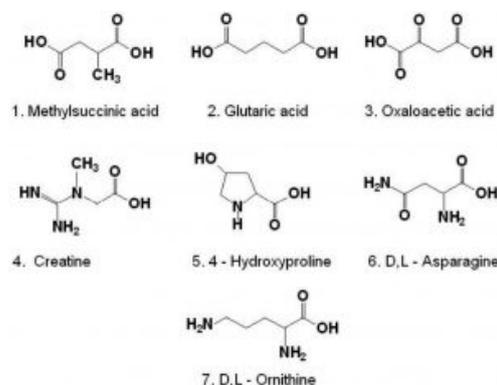


A Simple and Reproducible LCMS Method

This Method demonstrates the ability to separate a wide variety of polar metabolites using Mass Spectrometry for detection. When monitoring the EIC at 133 m/z, seven Peaks were positively detected in a Urine sample by comparison with standards or by using accurate mass.

Although primarily focused on the analysis of metabolites in Urine, the Method could also be applied to the determination of these compounds in other physiological or biological fluids.

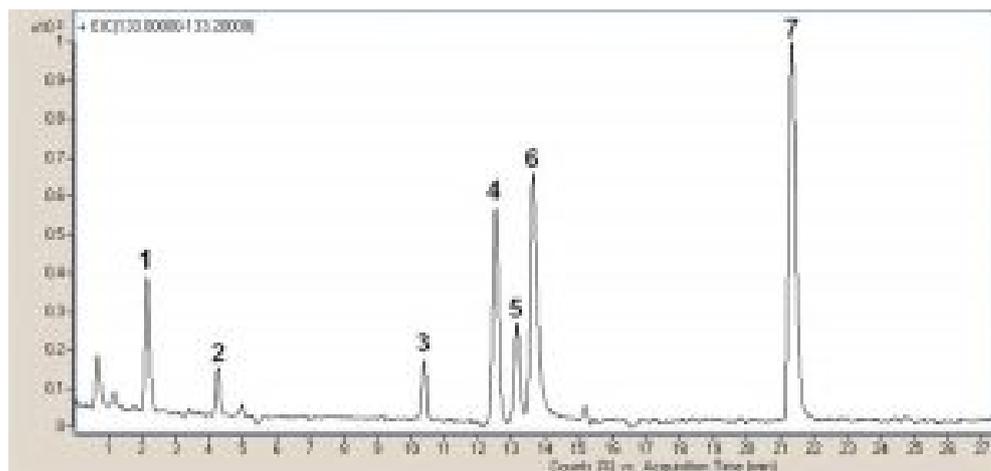

Peak:

1. Methylsuccinic Acid, 133 m/z (M+H)⁺, RT = 2.14 minutes
2. Glutaric Acid, 133 m/z (M+H)⁺, RT = 4.25 minutes
3. Oxaloacetic Acid, 133 m/z (M+H)⁺, RT = 10.40 minutes
4. Creatine, 133 m/z (M+H)⁺, RT = 12.52 minutes
5. 4-Hydroxyproline, 133 m/z (M+H)⁺, RT = 13.15 minutes
6. D,L-Asparagine, 133 m/z (M+H)⁺, RT = 13.62 minutes
7. D,L-Ornithine, 133 m/z (M+H)⁺, RT = 21.31 minutes

Method Conditions:

Column: Cogent Diamond Hydride™, 4μm, 100Å

Catalog No.: 70000-15P-2

Dimensions: 2.1 x 150mm

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Mobile Phase:

A: DI Water / 0.1% Formic Acid

B: Acetonitrile / 0.1% Formic Acid

Gradient:

Time (minutes)	%B
0	95
0.2	95
30	50
35	50
35.1	95

Post Time: 5 minutes**Flow Rate:** 0.4mL / minute**Detection:** ESI - POS - Agilent 6210 MSD TOF Mass Spectrometer

Sample Preparation: Synthetic Urine: 400µL of Acetonitrile was added to 100µL of Synthetic Urine and sample was centrifuged (3000g). Next 20µL of the supernatant was mixed with 10µL of 50:50 Acetonitrile / DI Water / 0.1% Formic Acid.

Notes: *The determination of metabolites and amino acids in biological fluids is an important problem in clinical biochemistry and analytical chemistry. Changes in the concentrations of these compounds in Urine, serum and other physiological fluids has proved to be correlated with several neurological disorders such as Alzheimer's disease, ischemic stroke as well as with a number of metabolic disorders such as phenylketonuria, argininemia, maple syrup Urine disease and others.*

Chromatogram presented is adapted from: J.J.Pesek, M.T. Matyska, S. M. Fisher, T. R. Sana, Journal of Chromatography A, 1204 (2008) p55.

**Attachment****No 88 Isobaric Metabolites in Urine Analyzed with HPLC pdf** 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.comWebsite: www.mtc-usa.com