


## Phosphorylated Sugars Analyzed with LCMS - AppNote

### Nucleotide Sugars UDP, ADP, TDP and CDP

UDP Hexanolamine was used as the Internal Standard in the analysis of Nucleotide Sugars in this rapid analysis Method. The Mobile Phase was designed to maximize the detector response in LCMS for maximum Chromatographic Efficiency. The simple “*Inverse Gradient*” which produces an Aqueous Normal Phase (ANP) HPLC Method, was required for the results shown.

This Method can be used in measuring Metabolite Concentration.

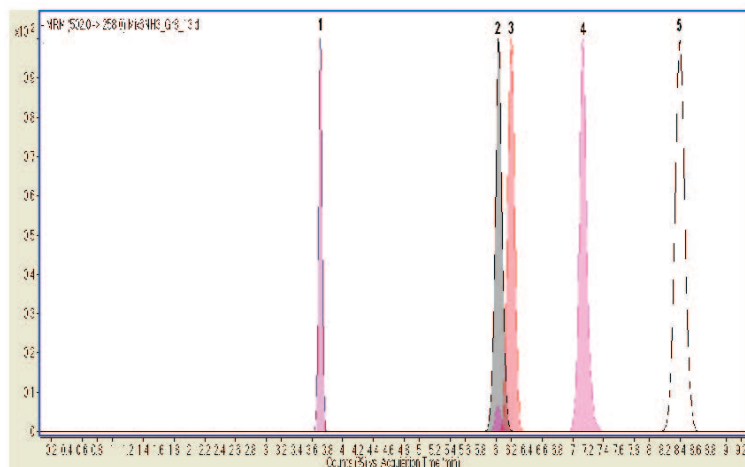



### Peaks:

1. ADP-Glucose, *RT* = 3.68 min, monitored MRM transitions were *m/z* 588 to *m/z* 346
2. Proprietary Sugar Nucleotide, *RT* = 6.03 min, monitored MRM transitions were *m/z* 563 to *m/z* 321
3. Proprietary Sugar Nucleotide, *RT* = 6.20 min, monitored MRM transitions were *m/z* 606 to *m/z* 385
4. CDP-Glucose, *RT* = 7.13 min, monitored MRM transitions were *m/z* 564 to *m/z* 322
5. UDP Hexanolamine (internal standard), *RT* = 8.40 min, monitored MRM transitions were *m/z* 502 to *m/z* 258

MRM = Multiple Reaction Monitoring in LCMS/MS

### Method Conditions

**Column:** Cogent Diamond Hydride™, 4µm, 100Å

**Catalog No.:** 70000-15P-2

**Dimensions:** 2.1 x 150mm

#### Mobile Phase:

A: DI Water / 0.1% Ammonium Formate *pH* 7.2

B: 90:10 Acetonitrile / DI Water / 0.1% Ammonium Formate *pH* 6

#### Gradient:

Time (minutes)	%B
0	95
10.0	75



12.0	75
12.1	95
15.0	95

**Flow rate:** 0.3mL / minute

**Sample Preparation:** 400mL DI Water / 400mL Acetonitrile / 20mL of Stock Solution of each compound / 5mL of 12% ammonia

**Detection:** ESI - NEG - Agilent 6210 MSD TOF Mass Spectrometer

**Note:** Please note the addition of small amount of Ammonia to the Sample Matrix. The alkaline environment of the Sample Matrix assured Efficient and Symmetrical Peaks for all analytes.



## Attachment

**No 70 Phosphorylated Sugars Analyzed with LCMS pdf** 0.3 Mb [Download File](#)