
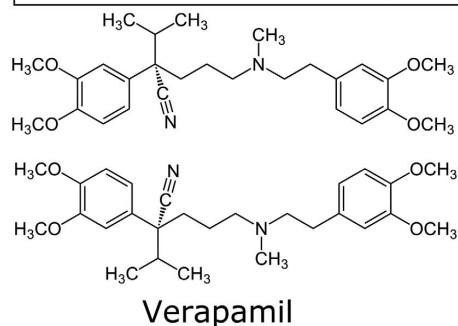


Alternative to USP Method with Improved Peak Efficiency

The USP impurities Method for Verapamil uses an Acetate Buffer with the ion pair agent 2-Aminoheptane to reduce Peak Tailing of the API and its Impurity. This Application Note uses an LCMS compatible Mobile Phase and produces Symmetrical Peaks for both Compounds. Excellent Resolution was also obtained between the two Peaks, far exceeding the USP System Suitability of $R_s \geq 1.5$.

This Method is a significant improvement over the USP monograph and illustrates the considerable potential for analyses of Amine containing solutes.

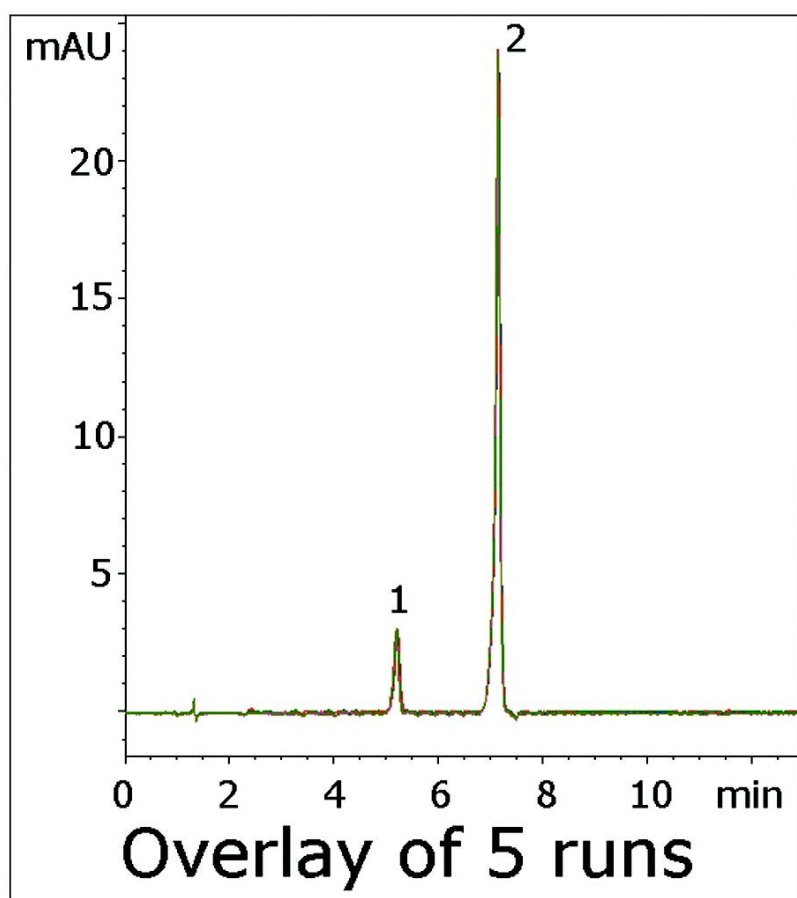

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

Peaks:

1. Verapamil Related Compound B
2. Verapamil

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Method Conditions

Column: Cogent Diamond Hydride™, 4µm, 100Å

Catalog No.: 70000-7.5P

Dimensions: 4.6 x 75mm

Mobile Phase:

A: DI Water / 10mM Ammonium Formate

B: 95:5 Acetonitrile / 10mM Ammonium Formate (v/v)

Gradient:

| Time (minutes) | %B |
|----------------|-----|
| 0 | 100 |
| 9 | 85 |
| 10 | 100 |

Post Time: 4 minutes

Injection vol.: 1µL

Flow rate: 1.0mL / minute

Detection: UV @ 278nm

Sample Preparation: Sample mixture contained 0.1mg / mL Verapamil and 0.02mg / mL Verapamil Related Compound B in a diluent of 50:50 Solvent A / Solvent B. Peak identities were confirmed with individual Standards.

t₀: 0.9 minutes

Note: Verapamil is an L-type calcium channel blocker used in the treatment of Hypertension, Angina Pectoris, and Cardiac Arrhythmia.



Attachment

No 180 Verapamil and Impurity Analyzed with HPLC pdf 0.5 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com