
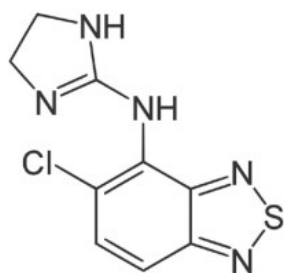


## Tizanidine HCl Tablet Analyzed by HPLC – AppNote

### Separation Method for Tizanidine Compatible with LC-MS

Tizanidine has numerous amine functional groups and can be a challenge for analysis by HPLC. The USP method uses Phosphate in the Mobile Phase which is not compatible with LC-MS. This Method however uses Formic Acid as the Mobile Phase additive and produces a sharp, symmetrical peak.

The USP system suitability for the tailing factor is not more than 1.6, and the Peak obtained has a value of 1.1. Data from two Column lots is shown in the figure, demonstrating Robustness of this Method.

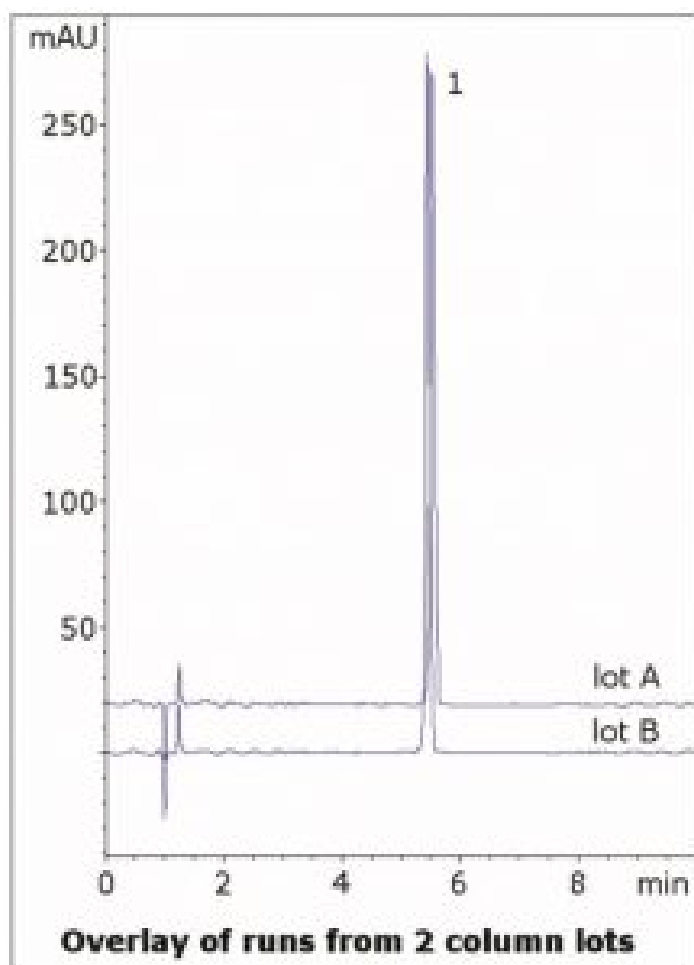

**Peak:**

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

## Method Conditions

**Column:** Cogent Diamond Hydride™, 4µm, 100Å

**Catalog No.:** 70000-7.5P

**Dimensions:** 4.6 x 75mm

**Mobile Phase:**

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

**Gradient:**

Time (minutes)	%B
0	95
1	95
6	40
7	95

**Post Time:** 3 minutes

**Injection vol.:** 1µL

**Flow rate:** 1.0mL / minute

**Detection:** UV @ 230nm

**Sample Preparation:** 4mg strength Tizanidine HCL tablet was ground and weighed in a 10mL volumetric flask. A portion of 50:50 Solvent A / Solvent B diluent was added and the flask was sonicated 10 minutes. It was then diluted to mark and filtered with a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.).

**t<sub>0</sub>:** 0.9 minutes

**Note:** Tizanidine is a centrally acting  $\alpha_2$ -adrenergic agonist used to treat spasms, cramping, tightness of muscles, and related conditions. It is available under the trade name Zanaflex® as well as generic versions.



### Attachment

No 232 Tizanidine HCL Tablet Analyzed by HPLC pdf 0.6 Mb [Download File](#)