


Compounds of the Carotenoid Family Retained

These compounds are very lipophilic and often have subtle differences in structure. Resolution of these compounds can be difficult with a standard Reversed Phase C8 or C18 Column because of these Similarities. The Cogent C30 Column on the other hand can further differentiate by analyte Shape in addition to Reversed Phase Interactions. At lower temperatures, the long alkyl chains become more rigid and steric effects become significant, leading to greater Selectivity.

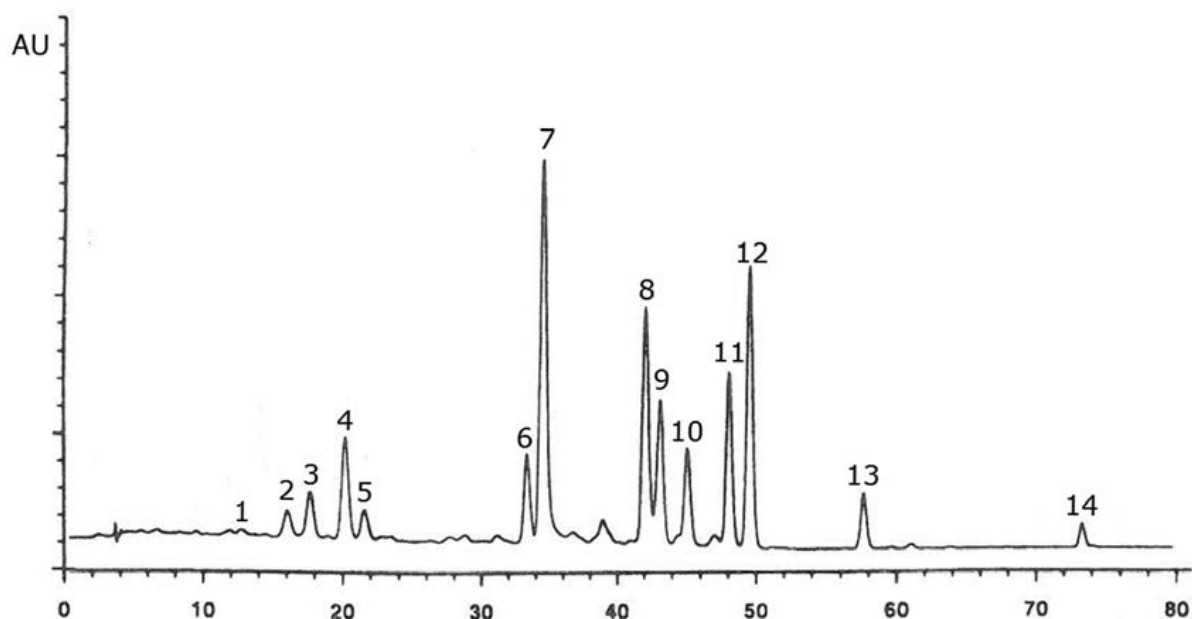

Peaks:

1. Astaxanthin
2. Capsanthin
3. Lutein
4. Zeaxanthin
5. Canthaxanthin
6. β -Kryptoxanthin
7. Echinenone
8. 15-cis β -Carotene
9. 13-cis β -Carotene
10. α -Carotene
11. trans β -Carotene
12. 9-cis β -Carotene
13. d-Carotene
14. Lycopene

Method Conditions:

Column: Cogent C30™, 5 μ m, 200Å

Catalog No.: 71030-25P

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Dimensions: 4.6 x 250mm

Mobile Phase:

A: 81:15:4 Methanol / Methyl Tert-Butyl Ether / DI Water (v/v)

B: 6:90:4 Methanol / Methyl Tert-Butyl Ether / DI Water (v/v)

Gradient:

Time (minutes)	%B
0	0
90	100

Temperature: 20°C

Flow rate: 1.0ml / minute

Detection: UV @ 450nm

Sample Preparation: Reference Standards of each Analyte in a mixture.

Notes: Carotenoids are a broad class of more than 600 compounds. They can be divided into two types: Xanthophylls and Carotenes. Xanthophylls contain oxygen in their structure and get their name from the Greek words Xanthos (yellow) and Phyllon (leaf). Carotenes are unsaturated hydrocarbons which do not contain other elements in their structure.



Attachment

No 262 Carotenoids Analyzed with HPLC pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com