
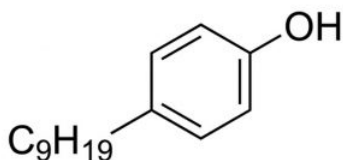


Nonylphenol Analyzed by HPLC - AppNote

Separation of Isomer Peaks

Nonylphenol is a very hydrophobic compound and is suitable for analysis by Normal Phase. It is produced commercially by Acid-catalyzed Alkylation of Phenol with a mixture of nonenes. Therefore a variety of product isomers are possible, with different branching of the C9 group and position of the chain on the ring.

Normal Phase HPLC is well-suited to isomer separations and in this method, separation was observed between the main peak and two smaller peaks. Three runs are shown to illustrate the repeatability, which is often a concern with Normal Phase Methods that use ordinary Silica Columns.

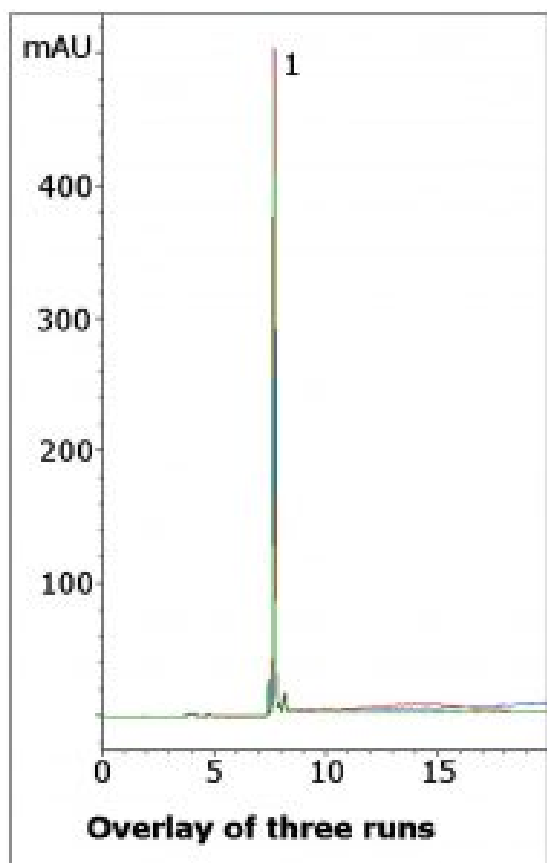

Nonylphenol

Main Peak:

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Method Conditions

Column: Cogent Silica-C™, 4µm, 100Å

Catalog No.: 40000-10P

MICROSOLV

Dimensions: 4.6 x 100mm

Mobile Phase:

A: Ethyl Acetate

B: Hexane

Gradient:

Time (minutes)	%B
0	100
4	100
19	90
20	100

Post Time: 3 minutes

Injection vol.: 1µL

Flow rate: 1.0mL / minute

Detection: UV @ 277nm

Sample Preparation: Nonylphenol reference standard dissolved in a Hexane diluent.

t₀: 1.3 minutes

Note: Nonylphenols are used in synthesis as a starting material for various surfactants. They are subjected to Ethoxylation to produce Alkyl-Phenol Ethoxylates.



Attachment

No 264 Nonylphenol HPLC Method pdf 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center
Copyright 2024, All Rights Apply
MicroSolv Technology Corporation
9158 Industrial Blvd. NE, Leland, NC 28451
tel. (732) 380-8900, fax (910) 769-9435
Email: customers@mtc-usa.com
Website: www.mtc-usa.com