


## Cytarabine assay by HPLC- AppNote

### A reliable assay for detection of a chemotherapeutic

Cytarabine belongs to the category of chemotherapy called antimetabolites. Antimetabolites are very similar to normal substances within the cell. When these substances are incorporated into the cellular metabolism, they are unable to divide. Antimetabolites are cell-cycle specific and target cells at very specific phases in the growth cycle.

This simple Isocratic method may be used to quantitate this important drug. The method maintains consistent and reliable results as the peak values were measured for week-long intervals. The %RSD for retention times for both standards and sample were less than 0.4%.

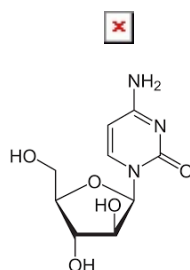

#### Peak:

Cytarabine

### Method Conditions:

**Column:** Cogent RP C18™ 5µm 100Å

**Catalog No.:** 68518-25P

**Dimensions:** 4.6 x 250mm

**Mobile Phase:** 5:95 methanol / buffer

*Buffer: 0.73g / L of monobasic sodium phosphate and 1.4g / L of dibasic sodium phosphate in DI water*

**Injection vol.:** 10µL

**Flow rate:** 1.0mL / minute

**Detection:** UV @ 254nm

**Sample Preparation:** Cytarabine 0.1mg / mL in DI water

**Note:** Cytarabine is an anti-cancer ("antineoplastic" or "cytotoxic") chemotherapy drug. This medication is classified as an "antimetabolite." Cytarabine is used to treat different forms of leukemia, some of these include acute and chronic myelogenous leukemia (AML and CML), acute lymphocytic leukemia (ALL), acute promyelocytic leukemia (APL). and also, Hodgkin's Lymphoma.

This method was developed by and is presented courtesy of **ARL- Eutech Scientific Services, Inc.**  
MicroSolv Technology Corporation



Printed from the Chrom Resource Center  
9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Date: 05-15-2024