
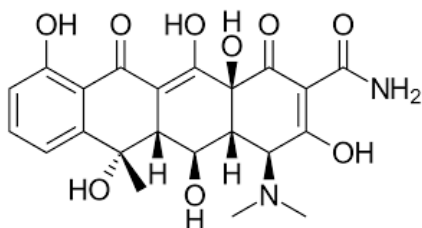


## Improved peak shape of a tetracycline antibiotic

Oxytetracycline has a role as an antibacterial drug, a protein synthesis inhibitor, an antimicrobial agent, an anti-inflammatory drug, and a bacterial metabolite. However, this compound interacts with stainless steel during HPLC and is known to have substandard peak shapes.

This data highlights important chromatographic performance improvements between our metal free treated stainless steel HPLC column and untreated stainless steel column with the same packing material and method conditions on the same instrument. Metal free surface treatments are available on all Cogent HPLC columns upon request.

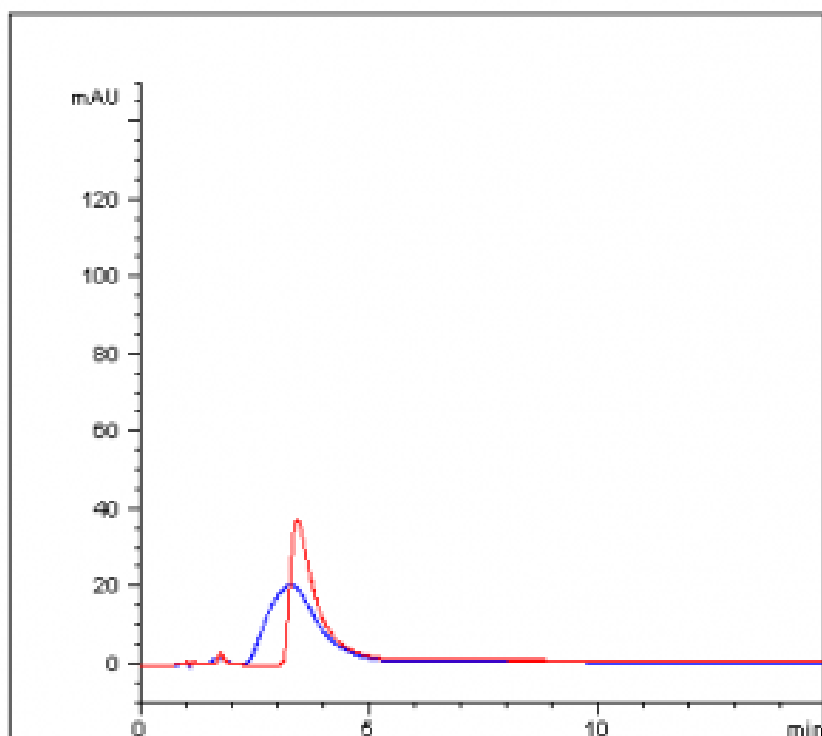

### Peak:

Oxytetracycline

Trace	Hardware	Width (Minutes)	Plates (N)
Red	Metal Free Coated Steel	0.46	304
Blue	non Coated Stainless Steel	1.18	43

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

**Method Conditions:****Column:** Cogent Diamond Hydride™, 4μm, 100Å**Catalog No.:** 70000-15D-2 and, 70000-15P-2**Dimensions of both columns:** 2.1mm x 150mm**Mobile Phase:** 50% acetonitrile / 50% DI water / 0.1% formic acid**Injection vol.:** 1μL**Flow rate:** 0.4mL / minute**Detection:** 270nm**Sample Preparation:** oxytetracycline 0.2mg / mL in DI water

*Note: Oxytetracycline is known as a broad-spectrum antibiotic due to its activity against such a wide range of infections. Oxytetracycline inhibits cell growth by inhibiting translation. It binds to the 30S ribosomal subunit and prevents the amino-acyl tRNA from binding to the A site of the ribosome.*



Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)Website: [www.mtc-usa.com](http://www.mtc-usa.com)