



Are the ends of the pre cut stainless steel tubing polished – FAQ


No, the ends of the pre-cut stainless-steel tubing are not polished. They are produced with a cutting device known as ECM (electro chemical machining) which produces a very flat, burr free edge which is better than polishing making them perfect for HPLC. See comparison of different types of cuts below.

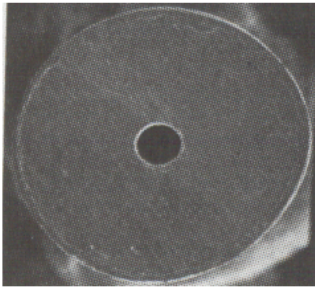

ECM-cut

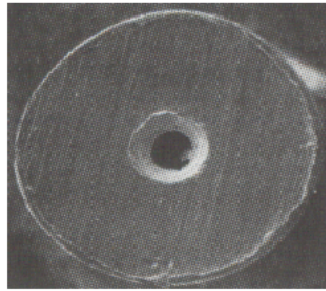

Mechanical cut

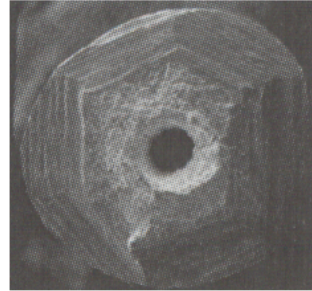

Manual cut

Click [HERE](#) for high pressure tubing ordering information and pictures.

Printed from the Chrom Resource Center

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Date: 05-19-2024