


Peptide Mixture Separated by HPLC - AppNote

Precise and Fast Equilibration Time

HPLC is a technique used for the characterization of Peptides and Reversed Phase HPLC is employed for the initial analysis and final, large scale purification. The first step of the production of Synthetic Peptides usually involves an initial Separation of these compounds from its mixture, on an analytical scale, then purification and collection of the target Peptide follows.

The Figure shown here presents the use of Reversed Phase HPLC Method in the Separation of a five Peptide Mixture. The 300Å pore size of the sorbent is ideal for separation of these Peptides chosen.

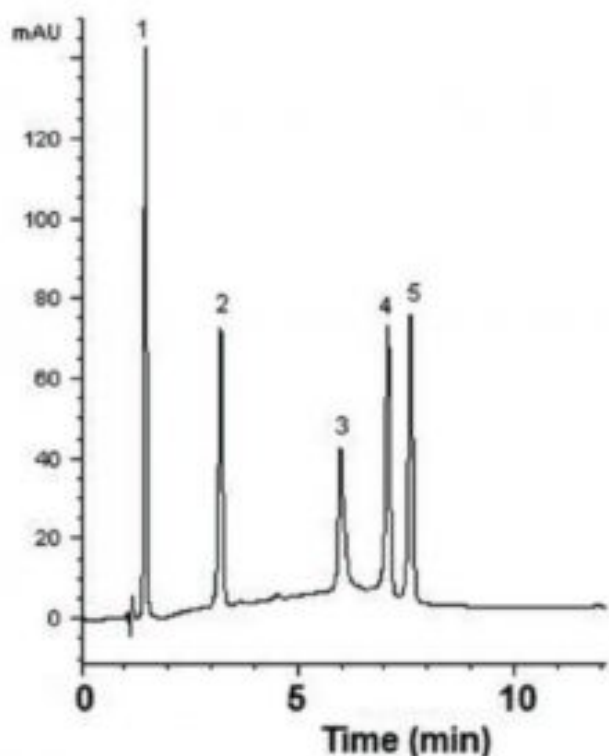

Peaks:

1. Gly-Tyr
2. Val-Tyr-Val
3. Methionine Enkephalin
4. Angiotensin II
5. Leucine Enkephalin

Method Conditions

Column: Cogent Bidentate C8 300™, 5µm, 300Å

Catalog No.: 40008-75P-3M

Printed from the Chrom Resource Center
Copyright 2024, All Rights Apply

MicroSolv Technology Corporation
9158 Industrial Blvd. NE, Leland, NC 28451
tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Dimensions: 4.6 x 75mm

Mobile Phase:

A: DI Water / 0.1% Trifluoroacetic Acid (TFA)

B: Acetonitrile / 0.1% Trifluoroacetic Acid (TFA)

Gradient:

Time (minutes)	%B
0	9
5	21
20	27
21	9

Post Time: 5 minutes

Flow rate: 1.0mL / minute

Detection: UV @ 214nm

Notes: Peptides are distinguished from proteins on the basis of the number of amino acid residues. Generally, this number is about 50 residues or fewer. Because of their smaller size, Peptides do not have the same degree of complexity that proteins do.



Attachment

No 91 Peptide Mixture Separation by HPLC pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com