



Furazolidone Analyzed with HPLC - AppNote

API and Related Compounds Separation

Furazolidone was Separated from related compounds in a commercially available pharmaceutical preparation. Excellent Peak shape and Selectivity from the related compounds was achieved with this Robust Method making it a great choice for QC applications even when using aggressive Mobile Phases that include acids. Molecular weight of Furazolidone is 225.16.

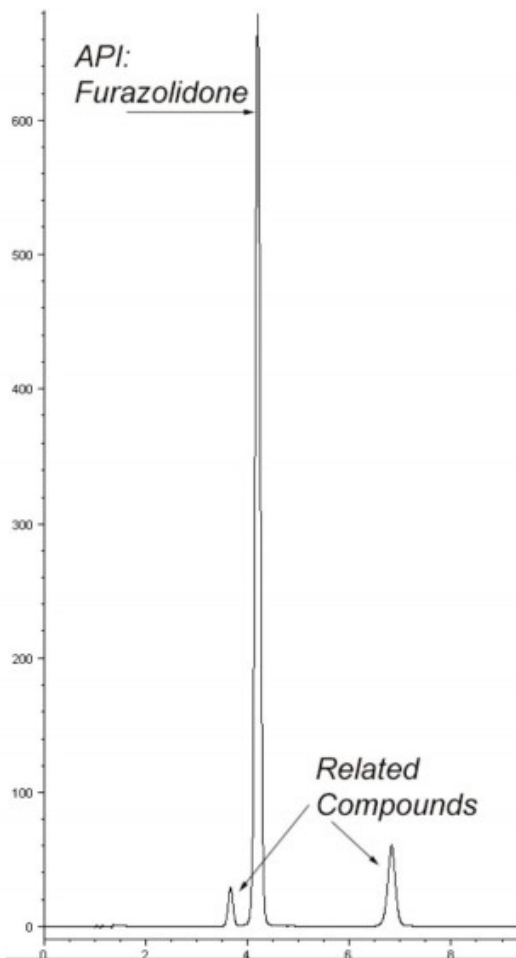

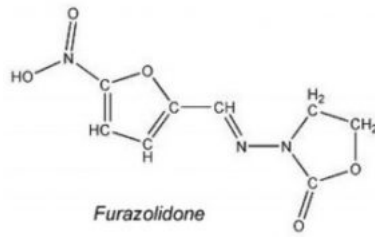

Peak:

Furazolidone

Method Conditions

Column: Cogent Bidentate C18™, 4µm, 100Å

Catalog No.: [40018-15P](#)

Dimensions: 4.6 x 150mm

Mobile Phase: 12% Acetonitrile / 88% DI Water / 1% Acetic Acid

Temperature: 30°C

Injection vol.: 20µL

Flow rate: 2mL / minute

Detection: UV @ 365nm

Note: This antibiotic drug is used to treat diarrhea and enteritis caused by bacteria or protozoan infection. Also used as a veterinary medicine for salmonids and *Myxobolus cerebralis* infection.



Attachment

No 73 Furazolidone Analyzed with HPLC pdf 0.2 Mb [Download File](#)