
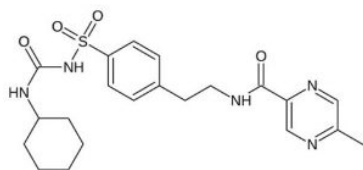


## Glipizide Analyzed with HPLC - AppNote

### Separation of API from Acid Degradation Products

This Method illustrates the capabilities to separate several Glipizide Degradants that are formed under acidic conditions. Good separation is obtained as well as Sharp, Symmetric Peak shapes for each compound.

The Retention times show good repeatability as well, as demonstrated in the Figure overlays. Also, the Mobile Phase Solvents are LCMS compatible, which expands the potential applications of the Method. Finally, Gradient Equilibration is fast and requires only 1 Column volume.

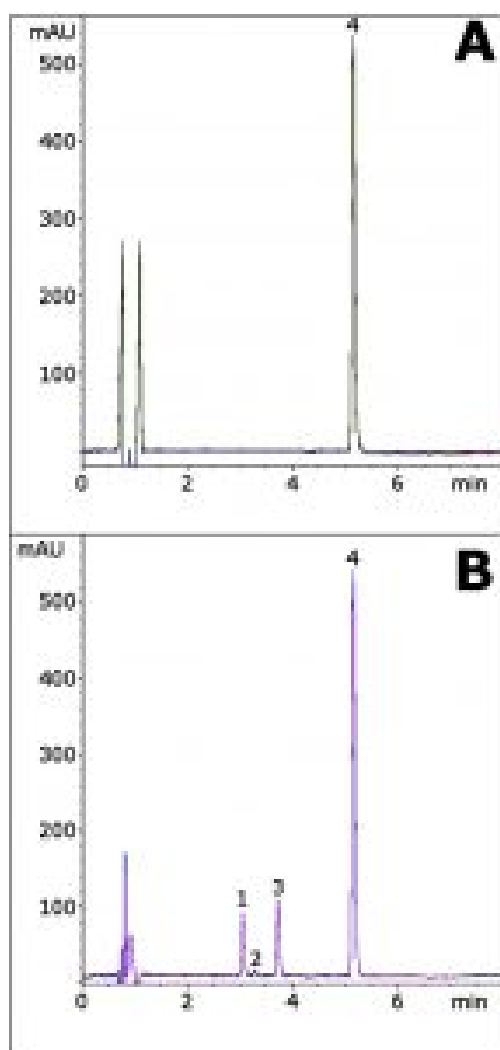

**Peaks:**  
1. Degradant

Printed from the Chrom Resource Center  
Copyright 2024, All Rights Apply  
**MicroSolv Technology Corporation**  
9158 Industrial Blvd. NE, Leland, NC 28451  
tel. (732) 380-8900, fax (910) 769-9435  
Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)  
Website: [www.mtc-usa.com](http://www.mtc-usa.com)

2. Degradant
3. Degradant
4. Glipizide

## Method Conditions

**Column:** Cogent Bidentate C18™, 4µm, 100Å

**Catalog No.:** 40018-75P

**Dimensions:** 4.6 x 75mm

**Mobile Phase:**

A: DI Water / 0.1% Formic Acid

B: Acetonitrile / 0.1% Formic Acid

**Injection vol.:** 10µL

**Flow rate:** 1mL / minute

**Detection:** UV @ 225nm

**Sample Preparation:**

*Figure. A (Non-Degraded):* 0.1mg / mL Glipizide in Methanol diluent.

*Figure. B (Acid Degradation):* 0.1mg / mL Glipizide in 50:50 Methanol /1N HCl diluent. Sample was heated at 85°C for 1 hour.

**t<sub>0</sub>:** 0.9 minutes

**Note:** Glipizide is an anti-diabetic sulfonylurea drug, sold under the trade name Glucotrol®.



## Attachment

**No 165 Forced Degradation of Glipizide with HPLC pdf** 0.3 Mb [Download File](#)