



## Coricidin HBP Day Tablet Analyzed by HPLC- AppNote

### Formulation of Guaifenesin and Dextromethorphan HBr


Click [HERE](#) for Column Ordering Information.

This Method shows the separation of two active ingredients in a formulation of Guaifenesin and Dextromethorphan. Dextromethorphan often exhibits tailing in many HPLC methods, but here an excellent Peak shape is observed.

This Method is suitable for routine assay of these compounds in combination formulations. A change in wavelength was used in order to obtain similar Peak heights between the two compounds, since their concentrations differ significantly.

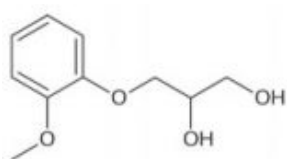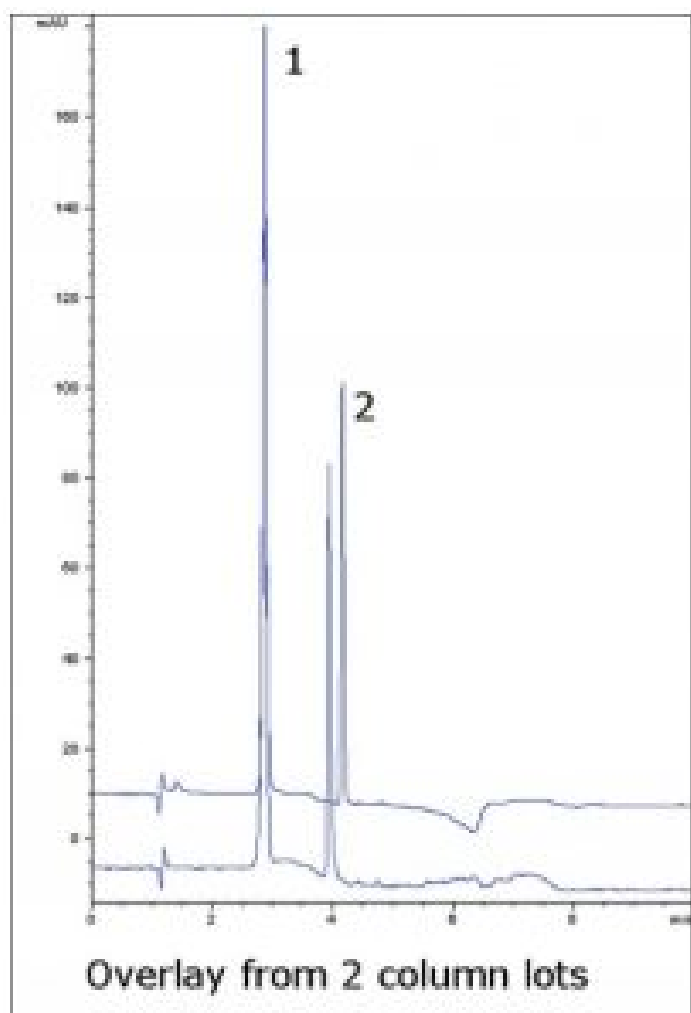

1

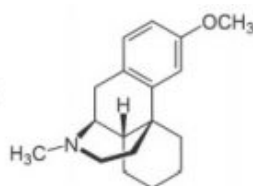

2

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

3 Industrial Blvd. NE, Leland, NC 28451

Phone: (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

1. Guaifenesin
2. Dextromethorphan

## Method Conditions

**Column:** Cogent Phenyl Hydride™, 4µm, 100Å

**Catalog No.:** 69020-75P

**Dimensions:** 4.6 x 75mm

**Mobile Phase:**

A: DI Water / 0.1% TFA (v/v)

B: Acetonitrile / 0.1% TFA (v/v)

**Gradient:**

Time (minutes)	%B
0	20
1	20
5	90
6	20

**Temperature:** 35°C

**Post Time:** 4 minutes

**Flow rate:** 1.0mL / minute

**Detection:** UV @ 254nm (0-3.5 minutes), 215nm (3.5-10 minutes)

**Injection vol.:** 2µL

**Sample Preparation:** Coricidin High Blood Pressure (HBP) Day® tablet containing 200mg Guaifenesin and 10mg Dextromethorphan HBr was ground and added to a 50mL volumetric flask. A portion of 50:50 Solvent A / Solvent B was added and the flask was sonicated 10 minutes. Then it was diluted to mark and mixed well. A portion was filtered with a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.). This filtrate was used for injections. The Peak identities were confirmed with individual standards.

**t<sub>0</sub>:** 0.9 minutes

**Note:** Guaifenesin is an expectorant agent and Dextromethorphan is an antitussive. This particular cold and cough formulation is designed for individuals with high blood pressure.



**Attachment**

**No 206 Coricidin HBP Day Tablet Analyzed by HPLC pdf** 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)