
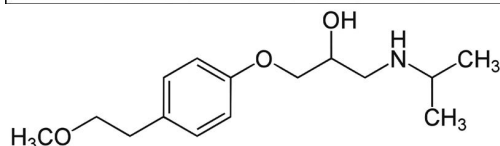


Metoprolol Tartrate Analyzed with HPLC – AppNote

Excellent Peak Shape for an Organic Amine

This Assay is easy to perform with a well-retained analyte Peak. Since Peak Tailing can often be an issue with organic amines, such as this compound, using conventional HPLC Columns it is notable that this Method produces a very Symmetrical Peak.

In the Chromatogram below, an overlay of injections from two lots of Columns is shown to illustrate Robustness and Precision of the Method which is suitable for routine assay of Metoprolol Tartrate Formulations.

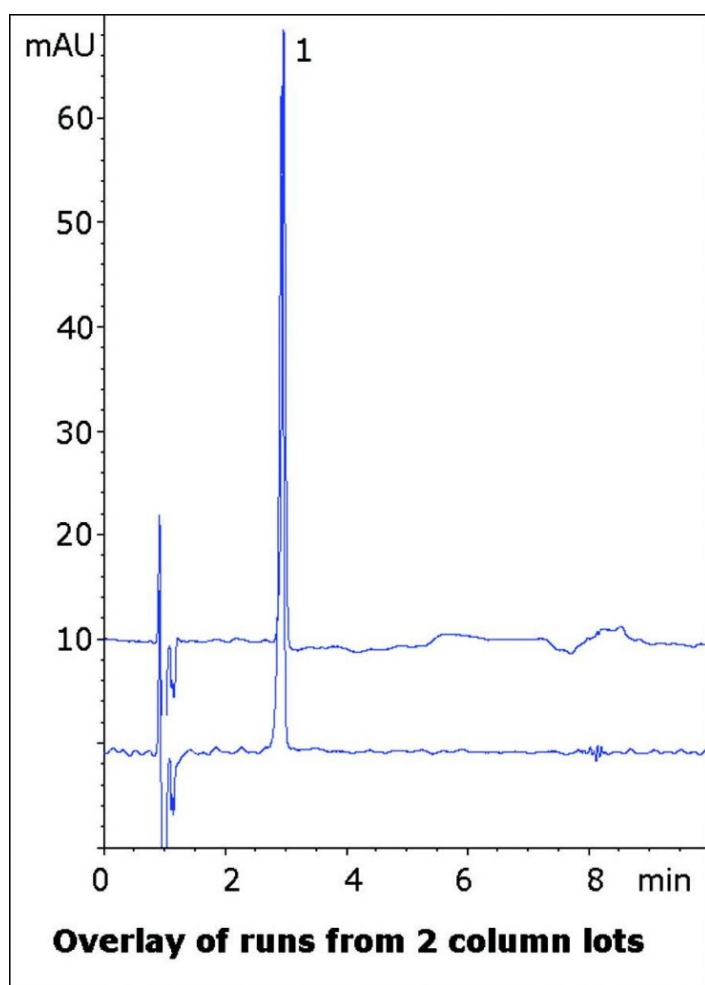

Peak:
Metoprolol

Method Conditions

Column: Cogent Diamond Hydride™, 4µm, 100Å

Printed from the Chrom Resource Center
Copyright 2024, All Rights Apply
MicroSolv Technology Corporation
9158 Industrial Blvd. NE, Leland, NC 28451
tel. (732) 380-8900, fax (910) 769-9435
Email: customers@mtc-usa.com
Website: www.mtc-usa.com

Catalog No.: 70000-7.5P

Dimensions: 4.6 x 75mm

Mobile Phase:

A: DI Water with 0.1% Trifluoroacetic Acid (TFA) v/v

B: Acetonitrile with 0.1% Trifluoroacetic Acid (TFA) v/v

Gradient:

Time (minutes)	%B
0	95
1	95
6	50
7	95

Post Time: 3 minutes

Injection vol.: 1µL

Flow rate: 1.0mL / minute

Detection: UV @ 215nm

Sample Preparation: 1mg Metoprolol Tartrate USP Reference Standard was dissolved in 1mL of 50:50 Solvent A / Solvent B. This stock solution was diluted 1:10 for HPLC injections using the same diluent.

t₀: 0.9 minutes

Note: Metoprolol is a selective β_1 receptor blocker. It is used to treat various cardiovascular conditions such as hypertension. Various trade names for the drug are available such as Lopressor® and Toprol®. It is sometimes found in combination formulations with hydrochlorothiazide.



Attachment

No 208 Metoprolol Tartrate Analyzed with HPLC pdf 0.4 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com