
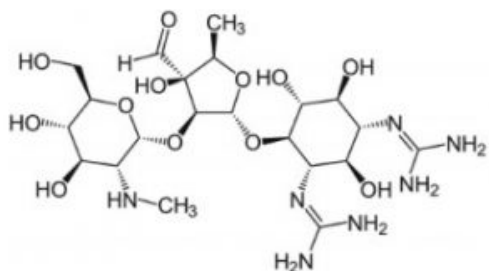


## Streptomycin Analyzed by HPLC – AppNote

### Retention of a Highly Polar Antibiotic

Streptomycin is very polar and difficult to retain with Reversed Phase methods using conventional Columns. In this Method, retention is achieved with Separation of the Streptomycin from two impurities found in the Standard Solution. Furthermore, the numerous amine groups do not cause tailing of Peak shape as with standard silica-based Columns.

Data from two Column lots shown in the figure illustrates the reproducibility of the Method and its robustness.

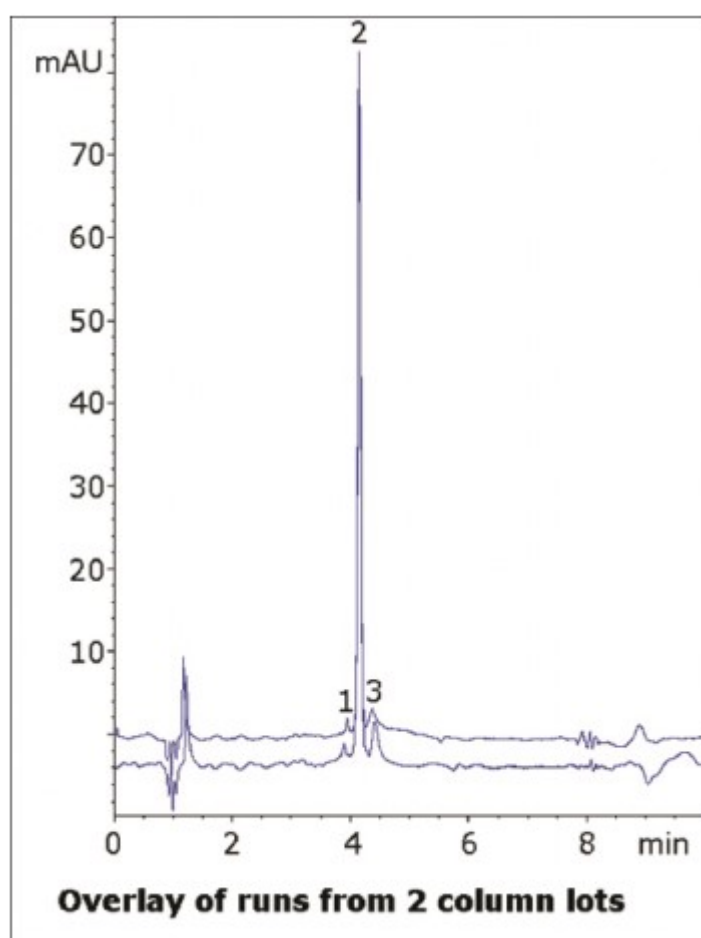

**Peaks:**  
1. Impurity

Printed from the Chrom Resource Center  
Copyright 2024, All Rights Apply  
**MicroSolv Technology Corporation**  
9158 Industrial Blvd. NE, Leland, NC 28451  
tel. (732) 380-8900, fax (910) 769-9435  
Email: customers@mtc-usa.com  
Website: www.mtc-usa.com

2. Streptomycin

3. Impurity

## Method Conditions

**Column:** Cogent Diamond Hydride™, 4µm, 100A

**Catalog No.:** 70000-7.5P

**Dimensions:** 4.6 x 75mm

### Mobile Phase:

A: DI Water / 0.1% Trifluoroacetic Acid (v/v)

B: Acetonitrile/ 0.1% Trifluoroacetic Acid (v/v)

### Gradient:

Time (minutes)	%B
0	95
1	95
6	40
7	95

**Post Time:** 3 minutes

**Injection Vol.:** 1µL

**Flow Rate:** 1.0mL / minute

**Detection:** UV @ 205nm

**Sample Preparation:** 1.0mg / mL Streptomycin Sulfate reference standard in Solvent A diluent.

**t<sub>0</sub>:** 0.9 minutes

**Note:** Streptomycin is an aminoglycoside antibiotic derived from the Actinobacterium *Streptomyces Griseus*. Streptomycin was isolated in 1943 and was the first antibiotic to treat tuberculosis.



Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)