
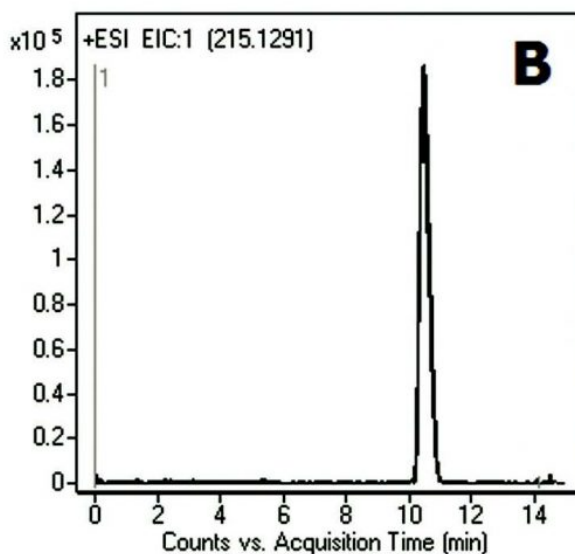


## Retention and Separation is Easy with this Method

3,3'-Diaminobenzidine (*DAB*) is a very challenging compound for analysis using HPLC with a UV detector. It is highly polar and hence difficult to retain when Reversed Phase HPLC Columns are used. Moreover, when there are a significant number of Silanol Groups present on the surface of the Column Stationary Phase, the Peak for *DAB* becomes very broad (5 - 10 minute Peak Width).

As can be seen from the Chromatograms below, this Method is an excellent choice for the analysis of *DAB* and resolves the issues of band broadening. The Peak shape is Symmetrical with high Efficiency.

The Repeatability of the Analysis is also remarkable as can be seen in a duplicate run shown in *Figure B*.

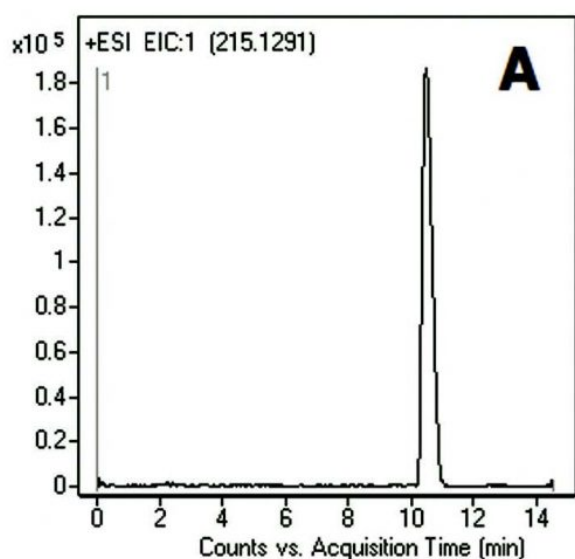

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply


**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

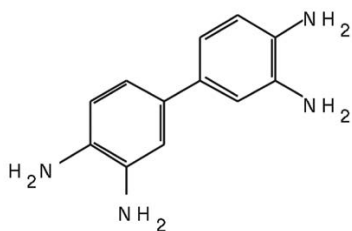**Peak:**3,3'-Diaminobenzidine 215.1291 m/z (M+H)<sup>+</sup>**Method Conditions:****Column:** Cogent Diamond Hydride™, 4µm, 100Å**Catalog No.:** 70000-15P-2**Dimensions:** 2.1 x 150mm**Mobile Phase:****A:** 50:50 DI Water / Methanol with 0.1% Formic Acid**B:** Acetonitrile with 0.1% Formic Acid**Gradient:**

Time (minutes)	%B
0	80
4	30
9	30
10	80

**Post Time:** 5 minutes**Injection vol.:** 1µL**Flow rate:** 0.4mL / minute**Detection:** ESI - POS - Agilent 6210 MSD TOF Mass Spectrometer**Sample preparation:***Working Solution:* Stock solution was diluted 1:100 with 50:50 Solvent A / Solvent B mixture.*Stock Solution:* 1mg / mL in DI Water Diluent. The Solution was filtered through a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.).**t<sub>0</sub>:** 0.9 minutes

**Notes:** DAB reacts with Hemoglobin (an oxidation reaction catalyzed by the Heme Groups) in the presence of Hydrogen Peroxide producing a dark brown color. This reaction is used to stain cells that were prepared with hydrogen peroxidase enzyme. DAB tablets are used in immuno-histology for the detection of Peroxidase sensitive Center Diaminobenzidine is a known mutagen (a compound that can induce changes in the genetic information of an organism).

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com



**Attachment**

**No 151 Diaminobenzidine Analyzed with LCMS pdf** 0.4 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)