
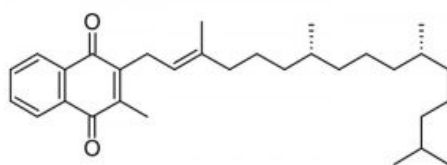


Vitamin K Isomers Analyzed by HPLC - AppNote

Phytonadione Separation by Shape Selectivity

In this Method, we separate the two E and Z isomers of Phytonadione on the basis of shape selectivity. The below chromatogram is a five injection overlay with a resolution value of 1.5.

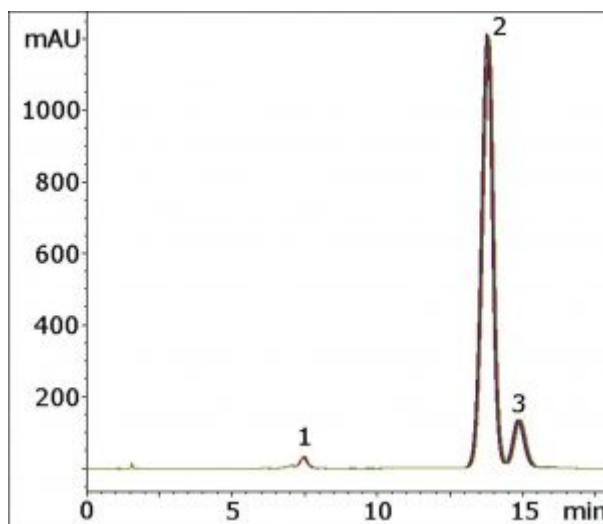

Peaks:

1. Impurity
2. Phytonadione (E isomer)
3. Phytonadione (Z isomer)

Method Conditions

Column: Cogent UDC-Cholesterol™, 4µm, 100Å

Catalog No.: 69069-15P

Dimensions: 4.6 x 150mm

Mobile Phase:

A: 50% DI Water / 50% MeOH / 0.1% Formic Acid

B: 97% Acetonitrile / 3% DI Water / 0.1% Formic Acid

Gradient:

| Time (minutes) | %B |
|----------------|----|
| 0 | 80 |
| 15 | 92 |
| 16 | 80 |

Temperature: 12°C

Post time: 2 minutes

Flow rate: 1.5mL / minute

Detection: UV @ 254nm

Sample Preparation:

Stock Solution: 10 μ L / mL Phytonadione in Acetonitrile diluent. (The solution was vortexed for 10 minutes.)

Working Solution: Stock solution was diluted 1:10 with Acetonitrile.

t₀: 1.0 minutes

Note: *Phytonadione (a.k.a. Phylloquinone, Vitamin K1) is a lipophilic vitamin that can be obtained in the diet from leafy green vegetables. It plays an essential role in blood clotting by acting as a cofactor for formation of coagulation factors II, VII, IX, and X. The letter designation for Vitamin K was based on the first letter of "Koagulationsvitamin" (coagulation vitamin), which is from the German journal that first published its identification by Danish biochemist Henrik Dam.*



Attachment

No 152 Vitamin K Isomers Analyzed by HPLC pdf 0.8 Mb [Download File](#)

Printed from the Chrom Resource Center

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Date: 07-22-2024