
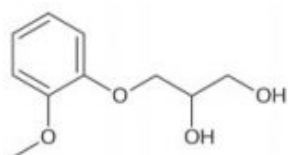


Formulation of Guaifenesin and Dextromethorphan HBr

Click [HERE](#) for Column Ordering Information.


This Method shows the separation of two active ingredients in a formulation of Guaifenesin and Dextromethorphan. Dextromethorphan often exhibits tailing in many HPLC methods, but here an excellent Peak shape is observed.

This Method is suitable for routine assay of these compounds in combination formulations. A change in wavelength was used in order to obtain similar Peak heights between the two compounds, since their concentrations differ significantly.

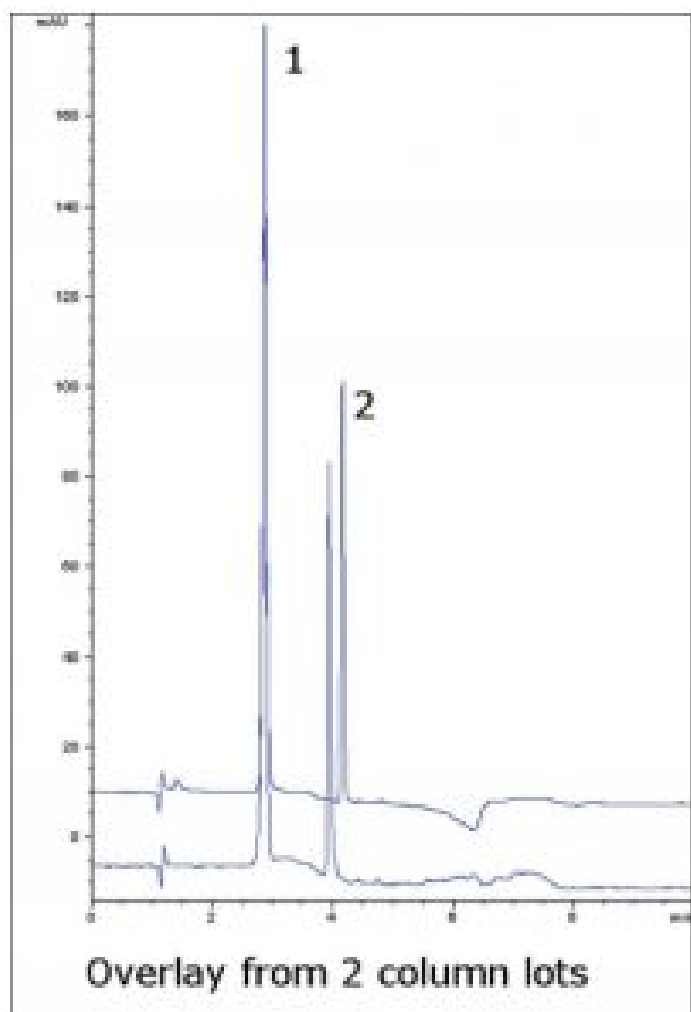

1

1. Guaifenesin
2. Dextromethorphan

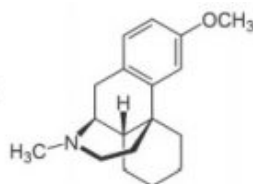

2

rinted from the Chrom Resource Center

Peaks: Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Method Conditions

Column: Cogent Phenyl Hydride™, 4µm, 100Å

Catalog No.: 69020-75P

Dimensions: 4.6 x 75mm

Mobile Phase:

A: DI Water / 0.1% TFA (v/v)

B: Acetonitrile / 0.1% TFA (v/v)

Gradient:

Time (minutes)	%B
0	20
1	20
5	90
6	20

Temperature: 35°C

Post Time: 4 minutes

Flow rate: 1.0mL / minute

Detection: UV @ 254nm (0-3.5 minutes), 215nm (3.5-10 minutes)

Injection vol.: 2µL

Sample Preparation: Coricidin High Blood Pressure (HBP) Day® tablet containing 200mg Guaifenesin and 10mg Dextromethorphan HBr was ground and added to a 50mL volumetric flask. A portion of 50:50 Solvent A / Solvent B was added and the flask was sonicated 10 minutes. Then it was diluted to mark and mixed well. A portion was filtered with a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.). This filtrate was used for injections. The Peak identities were confirmed with individual standards.

t₀: 0.9 minutes

Note: *Guaifenesin is an expectorant agent and Dextromethorphan is an antitussive. This particular cold and cough formulation is designed for individuals with high blood pressure.*



Attachment

No 206 Coricidin HBP Day Tablet Analyzed by HPLC pdf 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com