


## Your filter sparging stones with stems are listed as 1/8 or 1/16 ID shouldn't this be OD - FAQ

In the case of certain types of filter sparging stones where the tubing is inserted into the stone, the ID of the stems fits the OD of the low pressure tubing snugly.


The following schematic should clarify the above sentence. The ID of the sparging stone stem and the OD of the low pressure tubing are the same (**red** lines). The OD of the sparging stone stem (**green**) is larger than the ID of the low pressure tubing (**blue**).

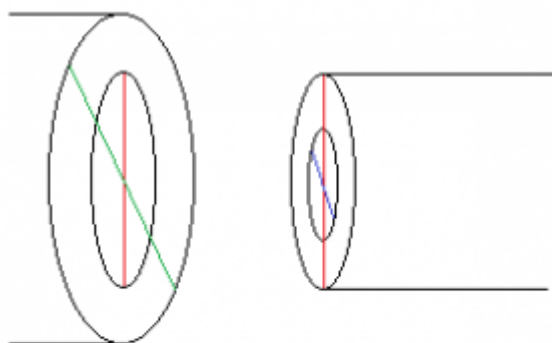

Sparging Stone Stem

Tubing

In another type of sparging stone, the OD of the sparging stone stem and the ID of the low pressure tubing are the same (**blue** lines). The ID of the sparging stone stem (**orange**) is smaller than the OD of the low pressure tubing (**red**).

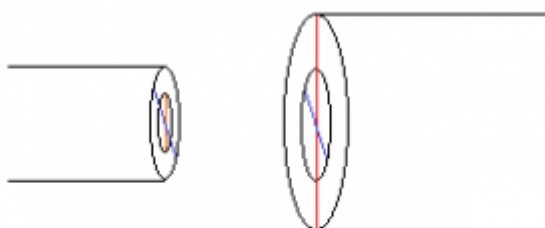

Sparging Stone Stem

Tubing

Click [HERE](#) for filter sparging stone ordering information and pictures.