
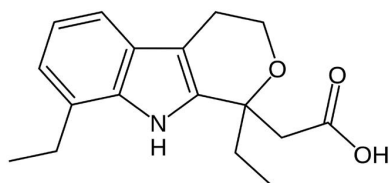


## Etodolac Tablet Analyzed with HPLC – AppNote

### An Assay Method for Etodolac Tablets

The USP Assay Method for Etodolac tablets uses a Phosphoric Acid based Mobile Phase which is not LCMS compatible. The Method in this Application Note is more versatile since it can be used for HPLC or LCMS. The Gradient can be adjusted if lower Retention is required.

Data from two Column lots is shown to illustrate Method Robustness and Consistency.

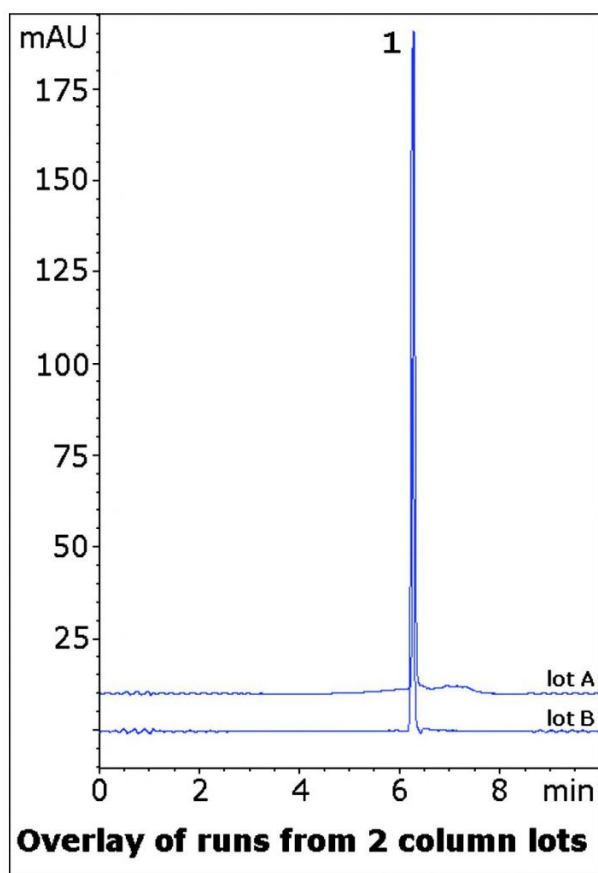

#### Peak:

Etodolac

### Method Conditions

**Column:** Cogent Bidentate C18™, 4µm, 100Å

**Catalog No.:** 40018-75P

**Dimensions:** 4.6 x 75mm

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

**Mobile Phase:**

A: DI Water with 0.1% Formic Acid (v/v)

B: Acetonitrile with 0.1% Formic Acid (v/v)

**Gradient:**

Time (minutes)	%B
0	10
1	10
6	80
7	10

**Post Time:** 3 minutes

**Injection vol.:** 1µL

**Flow rate:** 1.0mL / minute

**Detection:** UV @ 275nm

**Sample Preparation:** A portion of 50:50 Solvent A / Solvent B diluent was added and the flask was sonicated for 10 minutes. It was then diluted to mark and mixed. A portion was filtered with a 0.45µm Nylon Syringe Filter (MicroSolv Tech Corp.) and diluted 1:10.

**t<sub>0</sub>:** 0.9 minutes

**Note:** Etodolac is a prescription nonsteroidal anti-inflammatory drug. It is used for treatment of pain resulting from inflammation in conditions such as arthritis. It is sold in the USA under the trade name Lodine by Almirall Limited.

**Attachment**

No 247 Etodolac Tablet Analyzed by HPLC pdf 0.4 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)