

## Does MicroSolv have an equivalent vial to the Waters Total Recovery vial? – FAQ


The MicroSolv MRQ™ vials are designed for the same autosamplers as the Waters© Total Recovery Vial. However, this is not an exact match. The Total Recovery Vial is indicated for use when the HPLC or GC autosampler has an injection needle with a side port. The MicroSolv™ MRQ vials are made for when the autosampler has an injection needle with a bottom inlet or port.

The MRQ vials hold a little more liquid and have additional headspace. Both vials have the same style Well for maximum recovery of sample, both have the exact same 12 x 32 mm outer dimensions and screw threads, that are 9-425 mm. Both vials work very well with our caps.



Click [HERE](#) to view more details on the MRQ vials.

*Below is a “rough” photograph for comparison purposes only, of both the Total Recovery Vial (TRV) and an MRQ Vial. Note that the MRQ Vial has a larger total volume and a smaller bottom well. The TRV has a longer well and terminates much lower in the vial than the MRQ.*

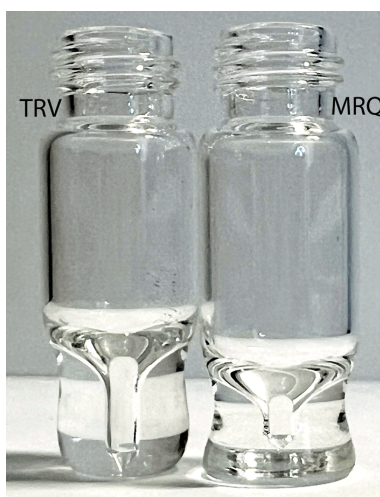

Printed from the Chrom Resource Center

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

Date: 07-27-2024