
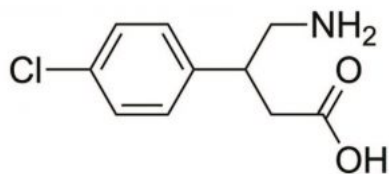


## Baclofen Analyzed with HPLC - AppNote

### Very Polar API Retained with Precision

This compound has two ionizable functional groups (an amine and a carboxylic acid), which makes it quite suitable for retention using this Method. A simple Gradient produced a very sharp Peak for this Analyte. These Method conditions can be applied to the routine Analysis of Baclofen Tablet formulations.

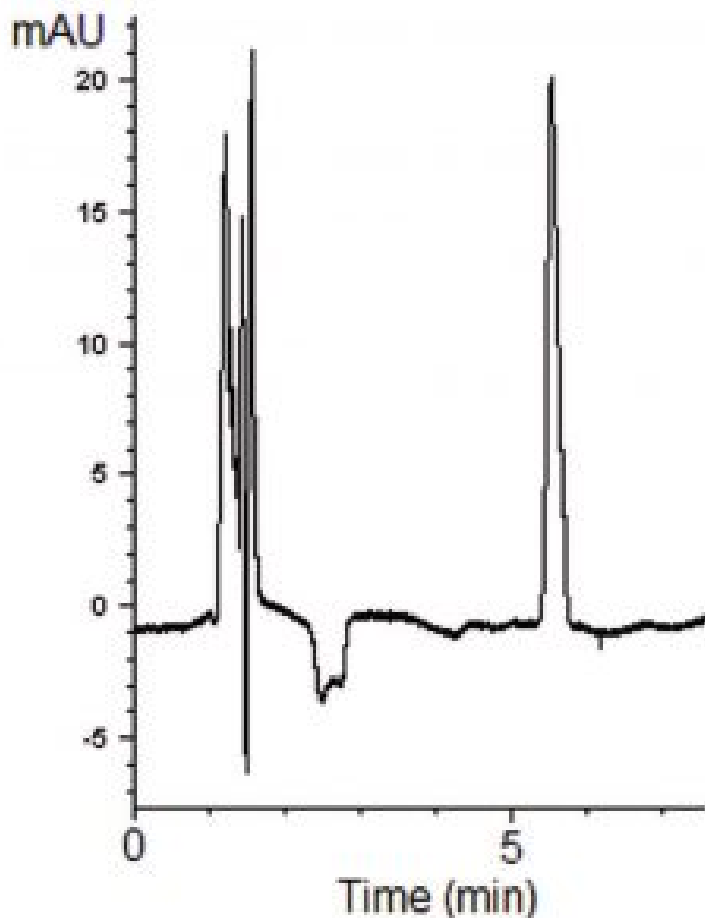

**Peak:**  
Baclofen

### Method Conditions

**Column:** Cogent Diamond Hydride™, 4µm, 100Å

**Catalog No.:** 70000-15P-2

**Dimensions:** 2.1 x 150mm

**Mobile Phase:**

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

A: DI Water with 0.1% Formic Acid (v/v)

B: Acetonitrile with 0.1% Formic Acid (v/v)

**Gradient:**

Time (minutes)	%B
0	95
5	30
6	30
7	95

**Injection vol.:** 1 $\mu$ L

**Flow rate:** 0.4mL / minute

**Detection:** UV @ 259nm

**Sample Preparation:** A Baclofen tablet (10mg strength) was ground with mortar and pestle, and 1.0mg of the ground powder was added to a beaker. Then 250 $\mu$ L diluent (66:4:30 DI Water / Glacial Acetic Acid / Methanol) was added and the solution was sonicated for 15 minutes. After gravity filtration, the filtrate was evaporated under vacuum and reconstituted in 250 $\mu$ L 90:10 Acetonitrile / DI Water with 0.1% Formic Acid (v/v).

*Note: Baclofen is used in the treatment of drug and alcohol addiction. In this usage, this drug lessens the symptoms of withdrawal. Treating addiction by prescribing baclofen is considered an "off label" use of this drug. It was originally designed to treat epilepsy. There is a need for an analytical method to monitor use of this medication.*



**Attachment**

**No 360 Baclofen Analyzed with HPLC pdf** 0.4 Mb

Printed from the Chrom Resource Center

Copyright 2024, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)