

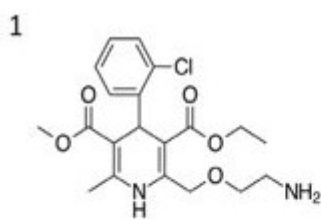


## Amlodipine, Valsartan, and Hydrochlorothiazide Tablet USP Assay - AppNote


This AppNote demonstrates the USP Method for Amlodipine, Valsartan, and Hydrochlorothiazide Tablets working with the Cogent RP C18. As shown in the Chromatograms, the Peak efficiency is superb and USP Peak tailing guidelines are easily met with the Column. This demonstrates a great alternative Column for easy plug and play with this USP Method.

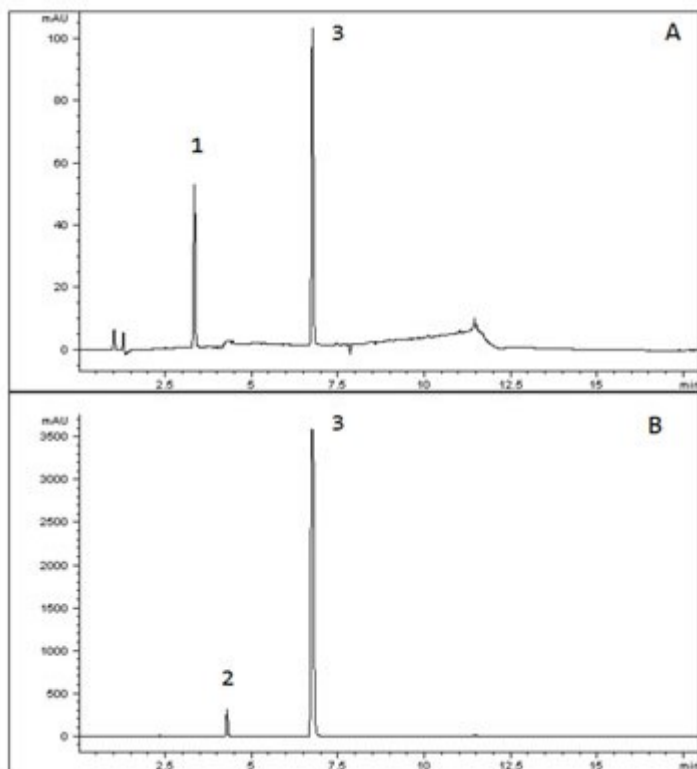

Amlodipine

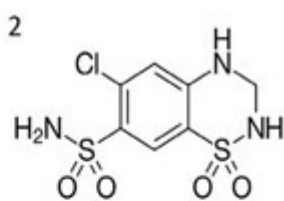

Hydrochlorothiazide

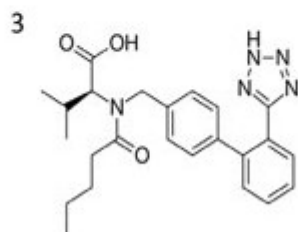

Valsartan

**Method Conditions:**

**Column:** Cogent RP C18™ 3um, 100A

**Catalog No.:** [68318-15P](#)

**Dimensions:** 4.6 x 150mm

**Mobile Phase:**

A. Acetonitrile/ DI Water / Phosphoric Acid (50:950:1)

B. Acetonitrile/ DI Water / Phosphoric Acid (950:50:1)

**Gradient:**

Time (minutes)	%B
0	5
3	50
6	60
10	95
10.1	5
15	5

**Injection Volume:** 10uL

**Flow Rate:** 1.5ml / minute

**Detection:** UV @ 225 nm

**Sample Preparation:**

A: Valsartan/Hydrochlorothiazide standards (0.16/0.025mg per mL)

B: Valsartan/Amlopidine Sample (3.2/0.1mg per mL)

USP Method Page Information: USP43-NF38 - 279 USP42-NF37 - 274 USP41-NF36 1S - 8278 Most Recently Appeared In: Pharmacopeial Forum: Volume No. 43(2)

**Note:** Amlodipine, Valsartan, and Hydrochlorothiazide is a combination of medicines that is used to treat high blood pressure (hypertension). Amlodipine is a calcium channel blocker (CCB) that regulates the absorption of calcium into the cells of the heart and blood vessels. This causes the blood vessels to relax and increases the supply of blood and oxygen to the heart while reducing its workload. Valsartan is an angiotensin II receptor blocker (ARB) and works by blocking a substance in the body that causes the blood vessels to tighten. This lowers blood pressure and increases the supply of blood and oxygen to the heart. Hydrochlorothiazide is a thiazide diuretic (water pill). It reduces the amount of water in the body.



**Attachment**

**N0.-381 Amlodipine, Valsartan, and Hydrochlorothiazide pdf** 0.1 Mb

Printed from the Chrom Resource Center

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

Date: 07-22-2024

[Download File](#)