
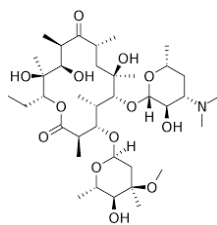


## Retention of Erythromycin analyzed with ELSD – AppNote

### Erythromycin, a macrolide antibiotic

As this compound lacks chromophores, it typically requires derivatization for use in UV detection. In this method, we retain excellent peak shape without the need for these pre-column derivatization steps saving time and resources.

RSD values (less than 0.5%) demonstrate the consistent and reliable retention, as shown in the 10 injections overlay below.

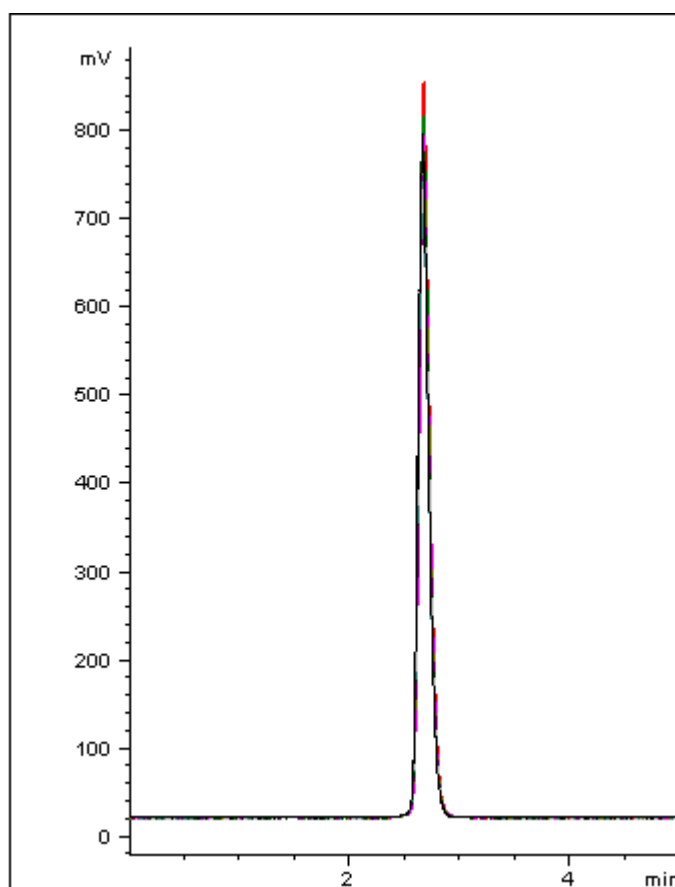

### Peak:

Erythromycin

### Method Conditions:

**Column:** Cogent Diamond Hydride™, 4μm, 100Å.

**Catalog No.:** [70000-10P](#)

# MICROSOLV

**Dimensions:** 4.6mm x 100mm

**Mobile Phase:** 50% acetonitrile / 50% DI water / 0.1% formic acid

**Injection vol.:** 1µL

**Flow rate:** 1.0mL / minute.

**Detection:** ELSD, gain: 9, temperature: 50°C, nitrogen: 3.5 bar.

**Sample Preparation:** 2.0mg / mL Erythromycin HCL in DI water.

*Notes: Erythromycin is a macrolide antibiotic used to treat bacterial infections and is also used to prevent recurrent rheumatic fever.*



Printed from the Chrom Resource Center

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)

Date: 05-15-2024